



Setting a new standard in von Willebrand Factor Activity Determination

INNOVANCE VWF Ac assay

[siemens-healthineers.us/hemostasis](https://www.siemens-healthineers.us/hemostasis)

SIEMENS
Healthineers



Von Willebrand disease: physiological causes and effects


Disturbances in the primary hemostasis system are the major cause of bleeding events. Patients with these disturbances are at increased risk for severe bleeding, especially during surgical interventions. Careful evaluation of primary hemostasis, which is principally determined by platelet function and von Willebrand factor activity, is important to minimize adverse outcomes in these patients.

Von Willebrand factor (VWF) is a multimeric, high-molecular glycoprotein involved in primary hemostasis, supporting platelet adhesion and aggregation under high-shear stress at the site of injury. Furthermore, VWF is the specific carrier protein of coagulation factor VIII, protecting factor VIII against inactivation and rapid clearance.

Von Willebrand disease (VWD), a dysfunction of primary hemostasis caused by quantitative or qualitative deficiency of von Willebrand factor, is the most common hereditary bleeding disorder. Leading symptoms of VWD include nosebleeds, menorrhagia, and bleeding after tooth extraction. In contrast, strong elevated VWF levels due to endothelial activation or dysfunction are frequently associated with cardiovascular disease and the risk for future ischemic heart disease or stroke.

Von Willebrand disease is present in roughly 1% of the general population. In patients undergoing surgery, a correct diagnosis allows targeted prophylactic treatment of the identified bleeding risk with the potential to significantly reduce the need for blood transfusions.¹


Leading symptoms of VWD

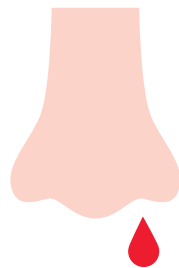

Nosebleeds



Bleeding after tooth extraction

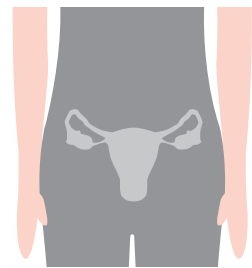

Menorrhagia

Distribution of VWF levels

A wide range of VWF levels is seen in the general population. There are several genes that affect VWF levels, with the major genetic influence coming from the ABO blood-group gene. Individuals with blood group O have 15–25% lower VWF levels than individuals with non-O blood group. However, a substantial fraction of the variation in VWF is not heritable and is related to inflammatory activity. VWF shows a fast and strong increase with acute phase reaction, which complicates the detection of mild VWD in the presence of stress or inflammation.²

NHLBI interpretation guide for VWF levels²

50–200% VWF	Normal VWF
30–50% VWF	Low VWF; probably modestly increased bleeding risk
<30% VWF	Decreased VWF; increased bleeding tendency; presence of VWF mutations likely

Subtypes of von Willebrand disease²

The most common type is VWF with quantitative defects, but a normal distribution of multimers.

VWD Type 1	Reduction of VWF defined; antigen and activity are in the same range due to decreased synthesis or increased degradation
VWD Type 3	Complete absence of VWF multimers (<5%)

Qualitative defects of VWF are less common.

VWD Type 2A	Selective reduction in or deficiency of large VWF multimers and reduction of platelet-dependent function
VWD Type 2B	Characterized by a strong affinity of VWF to GPIb, usually in conjunction with a low platelet count and reduced VWF activity
VWD Type 2M	Low platelet-dependent function of the VWF with large multimers present
VWD Type 2N	The FVIII binding site of VWF is defective; FVIII cannot bind and thus FVIII cannot be protected from premature degradation; clinical symptoms such as with hemophilia A

Management of VWD


To prevent or control bleeding in VWD patients, two major strategies are available. The first option is to increase plasma VWF concentration by stimulation and release of endogenous VWF stores in platelet and endothelial cells through application of desmopressin. This approach works well in most VWD type 1 patients and certain type 2 patients.

Replacement of VWF is the second alternative. Several human, plasma-derived VWF concentrates are available for therapy of those patients with insufficient response to desmopressin. However, before therapy can be started, a thorough anamnesis including appropriate lab testing is required.

Caveat: for VWD type 2B there is a clear contraindication.

Diagnosis and monitoring of von Willebrand disease

To correctly diagnose cases involving a bleeding anamnesis or a family history of von Willebrand disease and subsequently monitor patient response to therapy, physicians require a specific test for von Willebrand factor. A functional assay is the preferred screening method for detecting quantitative as well as qualitative defects of von Willebrand factor.³

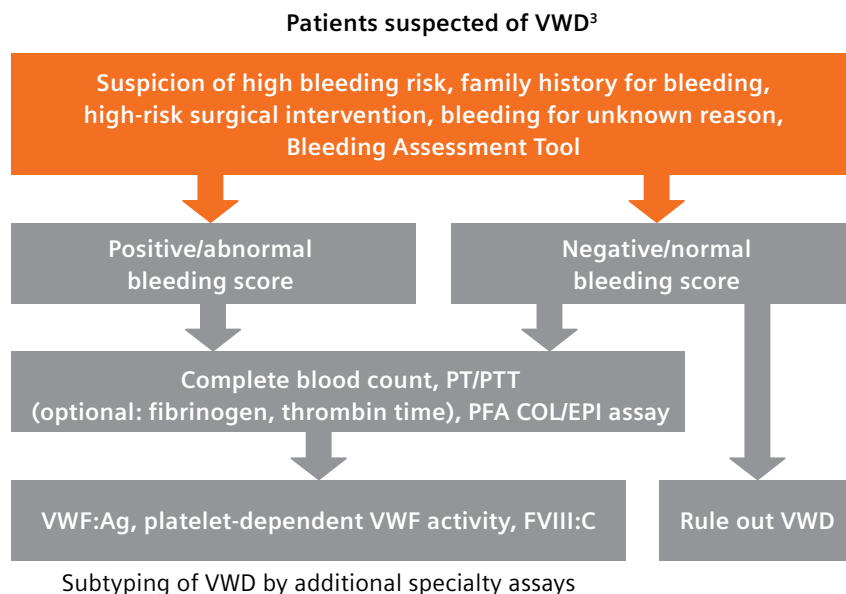

VWD subtype characteristics (modified from Reference 2)

Assay	Normal	Type 1	Type 2				Type 3
			2A	2B	2M	2N	
Screening tests							
PFA COL/EPI assay	N	N or ↑	↑	↑	↑	N	↑↑↑
Platelet count	N	N	N	N or ↓	N	N	N
VWF activity: • Platelet-dependent VWF activity	N	Borderline to ↓↓	↓ or ↓↓	↓↓	↓ or ↓↓	N or borderline	↓↓↓
Further VWF-related routine tests							
VWF Ag	N	Borderline to ↓↓	Borderline or ↓	Borderline or ↓	Borderline or ↓	N or borderline	↓↓↓
FVIII:C	N	N or ↓	N or ↓	N or ↓	N	↓↓	↓↓↓
VWF structure analysis							
VWF multimer analysis	Normal multimer distribution	Normal multimer distribution	Lack of HMWM	Lack of HMWM	Normal distribution	Normal multimer distribution	Absent
RIPA • High dose • Low dose	N NR	N or ↓ NR	↓ NR	N ↑	↓ NR	N NR	NR NR

N = Normal
NR = No reaction

A real functional assay for the automated determination of von Willebrand factor activity in human citrated plasma

INNOVANCE VWF Ac assay

The INNOVANCE® VWF Ac assay is a turbidimetric latex-based assay that requires no stirring and allows testing on all Siemens Healthineers and automated coagulation instruments.†

- Excellent precision and lot-to-lot consistency.
- Provides the ability to measure samples from 4–300% of norm based on three different settings in single-determination.
- Reagents are liquid and ready to use.

Assay performance	BCS® XP System	CS-2500 System	CS-5100 System
Required sample volume (µl)	45	16	16
Measuring range (% of norm)	4 - 300	4 - 300	4 - 300
Onboard stability (hours)†	36	24 / 72*	24 / 72*
Repeatability (CV %) of Control Plasma N	1.3	3	2.4
Within-Device/ Lab (CV %) of Control Plasma N	1.7	3.3	3

*Reagent I-II 24 hours, Reagent III 72 hours


INNOVANCE VWF Ac assay principle

Reagent I:
Latex particles coated with mouse anti-GPIb

Reagent II:
Buffer and blocking reagent

Reagent III:
Recombinant GPIb, carrying the gain-of-function mutations G233V and M239V

Sample:
Von Willebrand factor multimers

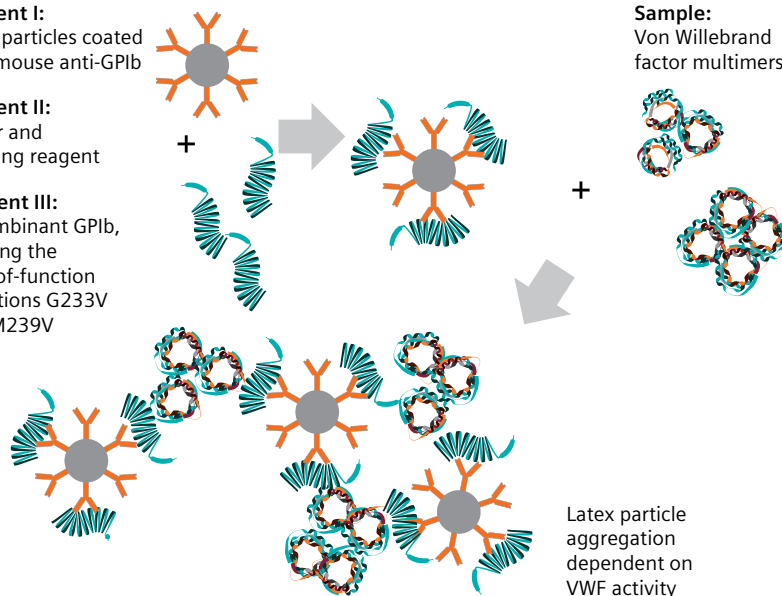

INNOVANCE VWF Ac assay employs an advanced new technology, allowing the assay to mimic the way in which VWF binds to glycoprotein Ib (GPIb), the major VWF receptor protein on platelets. Latex particles are coated with an antibody against GPIb, to which recombinant GPIb is added. The addition of patient plasma induces a VWF-dependent agglutination, which is detected turbidimetrically. Because the recombinant receptor protein includes two gain-of-function mutations, the assay does not require ristocetin.

†Except for CA-620/600 Systems.

Reliable, precise and sensitive

Bleeding disorders are a clinically severe complication requiring a reliable assay result to allow further therapeutic decisions within the shortest time.

High precision

The INNOVANCE VWF Ac assay provides excellent reproducibility over the full range of measurements.

	VWF Mean (% of norm)	Reproducibility (CV %)
Control Plasma N	87.75	2.95
Control Plasma P	28.97	5.30
Plasma pool (low)	9.14	3.48
Plasma pool (decision range)	49.26	2.83
Plasma pool (normal)	128.29	6.29
Plasma pool (high)	251.63	5.04

Multicentric reproducibility evaluation (5 days x 2 runs x 3 replicates design, one reagent/calibrator lot combination at three external sites with 89/90* samples).

*1 of 90 measurements (1.1%) was measured outside the measuring interval

Control and pool plasmas containing a range of 8–270% VWF were tested in two separate runs in two single determinations over 20 consecutive days without recalibration on all automated coagulation systems.

	Control N	Control P	Plasma pool (low)	Plasma pool (decision range)	Plasma pool (normal)	Plasma pool (high)
BCS XP System						
VWF Mean (% of norm)	84.83	28.62	10.06	52.23	129.74	271.47
Repeatability (CV %)	1.33	4.26	1.44	1.64	5.00	5.00
Within-Device/ Lab (CV %)	1.69	4.82	2.04	2.65	6.43	6.40
CS-2500 System						
VWF Mean (% of norm)	82.56	28.80	8.08	45.02	134.95	272.12
Repeatability (CV %)	3.03	2.38	2.20	2.14	3.44	2.60
Within-Device/ Lab (CV %)	3.33	2.63	2.41	3.52	4.25	3.06
CS-5100 System						
VWF Mean (% of norm)	81.46	30.16	8.82	46.44	133.98	253.08
Repeatability (CV %)	2.37	2.15	1.61	1.24	1.44	1.19
Within-Device/ Lab (CV %)	3.01	2.19	2.13	2.79	3.56	1.75


Diagnostic accuracy

Excellent correlation to VWF:RCo method

In the past, VWF:RCo assays such as BC von Willebrand Reagent were the recommended screening method for VWD. Ristocetin can also induce decapping of the A1 binding domain, the mechanism employed in VWF:RCo assays. However, ristocetin is an artificial substance, and in the presence of certain polymorphisms, false-low VWF levels may be measured that are not connected to bleeding risk.⁴ In contrast to VWF:RCo assays, assays such as the INNOVANCE VWF Ac assay directly determine real activity of VWF. Both assays show excellent correlation.

“The automated VWF:GPIbM has demonstrated excellent clinical performance and has been favored by experts as interchangeable with VWF:RCo.”

Higgins, et al. Am J Hematol. 2019;94:496-503.5

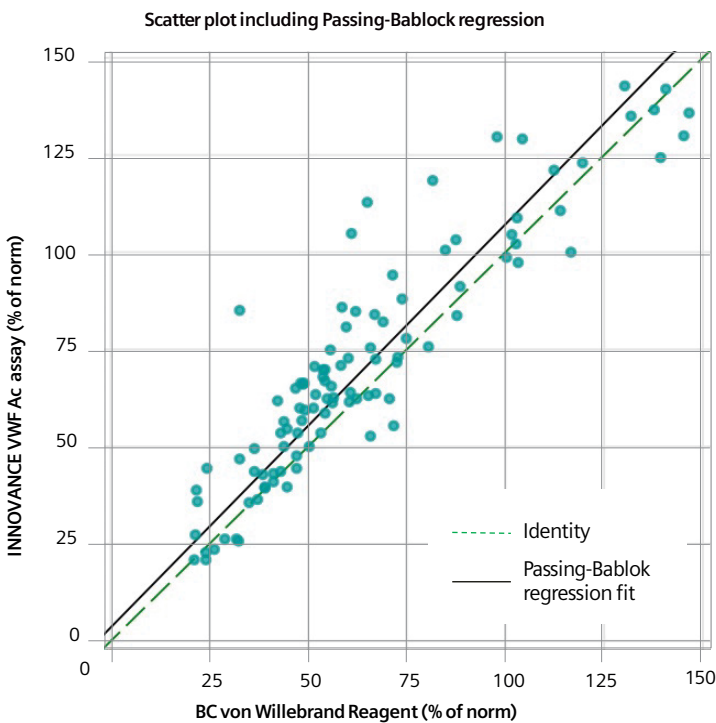

Method comparison measured on BCS® XP System

Number of samples (n)	102
Pearson Correlation Coefficient (r)	0.916
Slope	1.04
Intercept (% of norm)	3.50
Predicted Bias at 30% of norm	4.60% of norm
Predicted Bias at 50% of norm	10.14%

Siemens Healthineers offers a comprehensive portfolio of systems and assays for the diagnosis of bleeding disorders and therapy monitoring:

Platelet diseases

- PFA-100 System

Von Willebrand disease

- PFA-100 System
- BC von Willebrand Reagent
- INNOVANCE VWF Ac assay
- vWF Ag assay*

FVIII, IX deficiency, hemophilia A+B

- Coagulation Factor VIII Deficient Plasma
- Coagulation Factor IX Deficient Plasma
- Factor VIII Chromogenic Assay

Contact your Siemens Healthineers representative today for more information on our bleeding portfolio or explore more on:

➔ [siemens-healthineers.us/hemostasis](https://www.siemens-healthineers.us/hemostasis)

Diagnostic accuracy

Diagnostic accuracy of the INNOVANCE VWF Ac assay was evaluated in a method comparison study and a concordance study with patients representing the intended use population. Four clinical sites evaluated a total number of 144 patients with confirmed or suspected VWD. The patient cohort included 60 patients previously diagnosed with VWD [29 patients with VWD Type 1, 19 patients with VWD Type 2 (7 patients with Type 2A, 5 patients with Type 2B, 4 patients with Type 2M, 3 patients with Type 2N), 8 patients with Type 3, 4 patients with acquired VWD].

Additionally, 5 patients with haemophilia A, 1 patient with platelet dysfunction and 78 patients without final VWD related diagnosis at the time of enrollment were included.

After the exclusion of 6 patients due to therapy with increasing agents and diluted samples, an independent physician interpreted the laboratory values and classified the remaining patients (n=138) patients into one of the following categories: VWD, Low VWF, VWD Excluded, or Inconclusive Diagnosis.

Concordance Evaluation		Laboratory classification using the Standard of Care (SOC) VWF activity device in the initial VWD test panel				
		VWD (either type)	"low" VWF acc. to NIH guideline	VWD excluded	Inconclusive diagnosis	Total
Laboratory classification using the INNOVANCE VWF Ac assay in the initial VWD test panel	VWD (either type)	26	1	0	0	27
	"low" VWF acc. to NIH guideline	3	10	2	0	15
	VWD excluded	2	7	68	2	79
	Inconclusive diagnosis	4	1	1	11	17
	Total	35	19	71	13	138

Patients for whom the classifications of the initial VWD test panel results were discordant were further assessed: The discordances were evaluated in a Second Line Assessment. This assessment included an evaluation of potential benefits or risks of using either the SOC VWF activity device or the INNOVANCE VWF Ac assay in the initial VWD test panel.

VWF activity assay used in the initial test panel	No Impact	Potential Benefit	Potential Risk	Potential Benefit and Potential Risk
INNOVANCE VWF Ac	12	9	3	1*
SOC	14	2	8	1*

*The diagnostic classification/interpretation of one was evaluated as potentially beneficial and potentially risky for both VWF activity methods.

Positive benefit-risk-analysis

Results

The results demonstrate that the benefit of using INNOVANCE VWF Ac assay is higher compared to using the SOC VWF activity assay whereas the risk of using INNOVANCE VWF Ac assay is lower compared to the SOC VWF activity assay.

Conclusion


3-fold higher number of benefits associated with the use of INNOVANCE VWF Ac assay!

Be state of the art

The ASH/ISTH/NHF/WFH 2021 guideline on the diagnosis of von Willebrand disease suggests the use of VWF:GPIbM and VWF:GPIbR assay over the VWF:RCo assay



'The panel suggests newer assays that measure the platelet binding activity of VWF (eg, VWF:GPIbM, VWF:GPIbR) over the VWF:RCo assay (automated or non-automated) for the diagnosis of VWD [...].'

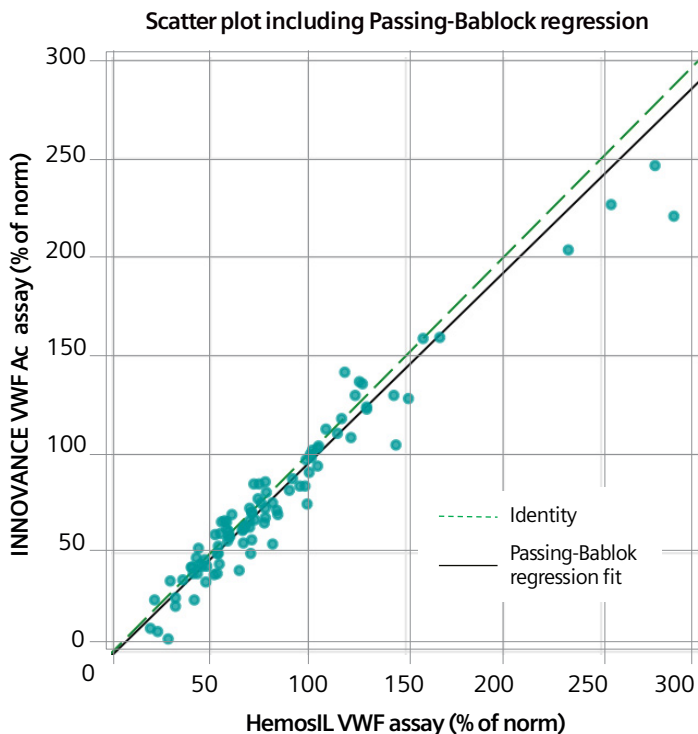

Good comparison to automated competitive assays


In the U.S., there are only a few automated von Willebrand activity assays available.

The HemosIL von Willebrand Factor Activity assay is a latex particle enhanced immunoturbidimetric assay based on the binding of a specific anti-VWF antibody against the platelet binding site of VWF.

Number of samples: 97

Pearson correlation coefficient: 0.972
 $y=0.969x - 1.495$

"After the positive experience we had using this assay in a study, we as caregivers in the U.S. are urgently waiting to use the INNOVANCE VWF Ac assay to diagnose our patients as we believe it supports improved von Willebrand disease (VWD) diagnosis and patient management, particularly now given the new recommendations from the ASH ISTH NHF WFH 2021 guidelines on the diagnosis of VWD."

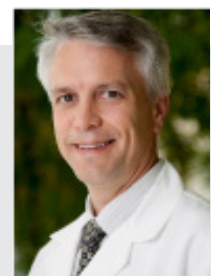

Steven Pipe, MD,
Professor of Pediatrics and Pathology
Laurence A. Boxer Research Professor of Pediatrics and Communicable Diseases
Pediatric Medical Director, Hemophilia and Coagulation Disorders Program
Director, Special Coagulation Laboratory
University of Michigan, MI US

Convenient and economical

Bleeding patients are always emergency cases who frequently come in at night or during the weekend.

Fast test results are required to support timely therapeutic decisions. The INNOVANCE VWF Ac assay is an automated, easy-to-use method available 24 hours a day/7 days a week to help physicians make the correct diagnosis quickly.


Convenient and easy to use

The test features liquid, ready-to-use reagents, with no processing preparation steps required for fully automated testing.

Stable calibration and control recovery

Controls Plasma N and P as well as Control Plasma P 1:3 diluted were tested over approx. 150 days without recalibration; all control results were well within the target range. Data from internal study.

Stable calibration of 8 month or longer with excellent control recovery minimizes recalibration requirements and increases availability of VWD testing.

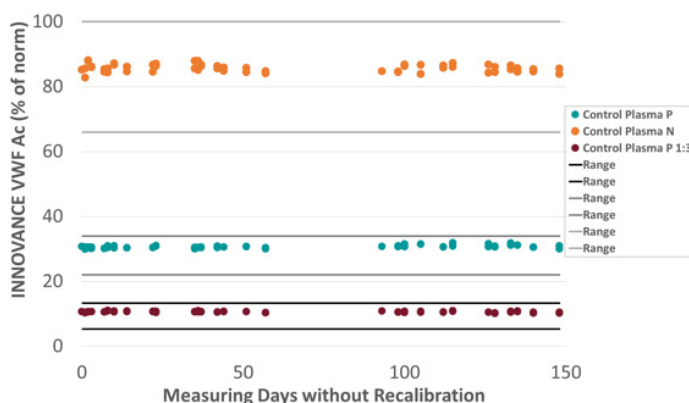

Stable reagents

Excellent reagent stability minimizes waste and maximizes the number of tests obtainable from a given volume of reagent:

Stability at 2–8°C once opened (BCS XP System)

Reagent I and II	37 days
Reagent III	113 days

Onboard stability

	BCS XP	CS-Systems
Reagent I and II	36 hours	24 hours
Reagent III	36 hours	72 hours

If the reagents are removed and stored closed at 2 to 8 °C in between the measurement periods, the maximum accumulated onboard stability time on rack line 1 to 4 on the BCS XP System is 24 hours for INNOVANCE VWF Ac REAGENT I and INNOVANCE VWF Ac REAGENT II and 36 hours for INNOVANCE VWF Ac REAGENT III.

(see Instruction for Use or Application sheet for further information)

Stable, robust method with low influence from interfering substances

With the INNOVANCE VWF Ac assay, no interference is observed up to:

- 656% VWF with regard to high-dose hook effect
- 1000 mg/dL hemoglobin
- 60 mg/dL bilirubin (unconjugated)
- 40 mg/dL bilirubin (conjugated)
- 726 mg/dL Lipids evaluated with lipemic samples
- 438 IU/mL rheumatoid factor

Available for labs of all sizes

The INNOVANCE VWF Ac assay makes it easy to perform comprehensive, highly accurate VWF testing.

Availability of the test across different platforms makes around-the-clock VWF Ac testing possible for labs of all sizes, eliminating the need to outsource testing.

The assay is available for:

- BCS XP Systems
- CS-2500 System
- CS-5100 System

Economical to use

The INNOVANCE VWF Ac assay is economical to use because of its:

- Single-run determination with excellent precision
- Low recalibration frequency (up to 8 month calibration stability)
- Long reagent stability
- Long shelf-life (24 months)



Siemens Healthineers AG (listed in Frankfurt, Germany: SHL) pioneers breakthroughs in healthcare. For everyone. Everywhere. Sustainably. As a leading medical technology company headquartered in Erlangen, Germany, Siemens Healthineers and its regional companies are continuously developing their product and service portfolio, with AI-supported applications and digital offerings that play an increasingly important role in the next generation of medical technology. These new applications will enhance the company's foundation in in-vitro diagnostics, image-guided therapy, in-vivo diagnostics, and innovative cancer care.

Siemens Healthineers also provides a range of services and solutions to enhance healthcare providers' ability to provide high-quality, efficient care. In fiscal 2022, which ended on September 30, 2022, Siemens Healthineers, which has approximately 69,500 employees worldwide, generated revenue of around €21.7 billion and adjusted EBIT of almost €3.7 billion.

Further information is available at www.siemens-healthineers.com.

The outcomes and statements provided by customers of Siemens Healthineers are unique to each customer's setting. Since there is no "typical" hospital and many variables exist (e.g., hospital size, case mix, and level of service/technology adoption), there can be no guarantee that others will achieve the same results.

On account of certain regional limitations of sales rights and service availability, we cannot guarantee that all products included in this brochure are available through the Siemens Healthineers sales organization worldwide. Availability and packaging may vary by country and is subject to change without prior notice. Some/All of the features and products described herein may not be available in the United States.

The information in this document contains general technical descriptions of specifications and options as well as standard and optional features, which do not always have to be present in individual cases.

Siemens Healthineers reserves the right to modify the design, packaging, specifications, and options described herein without prior notice. For the most current information, please contact your local sales representative from Siemens Healthineers.

Note: Any technical data contained in this document may vary within defined tolerances. Original images always lose a certain amount of detail when reproduced.

BCS, INNOVANCE, PFA-100, and all associated marks are trademarks of Siemens Healthcare Diagnostics Inc., or its affiliates. All other trademarks and brands are the property of their respective owners.

References:

1. Koscielny J, von Tempelhoff GF, Ziemer S, et al. A practical concept for preoperative management of patients with impaired primary hemostasis. *Clin Appl Thrombosis/Haemostasis*. 2004;10:155–66.
2. Nichols WL, Hultin MB, James AH, et al. von Willebrand disease (VWD): evidence-based diagnosis and management guideline, the National Heart, Lung, and Blood Institute (NHLBI) Expert Panel report (USA). *Haemophilia*. 2008;14:171–232.
3. James PD, et al. ASH ISTH NHF WFH 2021 guidelines on the diagnosis of von Willebrand disease. 2021 Jan 12 [published online]. DOI 10.1182/
4. Flood VH, Gill JC, Morateck PA, et al. Common VWF exon 28 polymorphisms in African Americans affecting the VWF activity assay by ristocetin cofactor. *Blood*. 2010;116:280–6. [bloodadvances.2020003265](https://doi.org/10.1182/blood-2009-08-20003265)
5. Higgins, et al. Automated assays for von Willebrand factor activity. *Am J Hematol*. 2019;94:496–503.

* Pending 510(k) clearance. Not available for sale in the U.S.

Siemens Healthineers Headquarters

Siemens Healthineers AG
Siemensstr. 3
91301 Forchheim, Germany
siemens-healthineers.com

USA

Siemens Healthcare Diagnostics Inc.
Laboratory Diagnostics
511 Benedict Avenue
Tarrytown, NY 10591-5005, USA
siemens-healthineers.us