

**White paper**

# **Mean Glandular Dose with MAMMOMAT Fusion**


[siemens-healthineers.com/mammomat-fusion](https://siemens-healthineers.com/mammomat-fusion)



## Introduction

MAMMOMAT Fusion\* is a premium mammography system made to enhance everyday screening and diagnostics. This includes intelligent dose reduction features and a new generation Cesium-Iodide (CsI) high-resolution detector to reliably deliver excellent images at the right dose – especially important when screening healthy women.

This paper presents the mean glandular dose with MAMMOMAT Fusion compared to a previous study<sup>1</sup> with other Full Field Digital Mammography (FFDM) systems in the clinical routine.

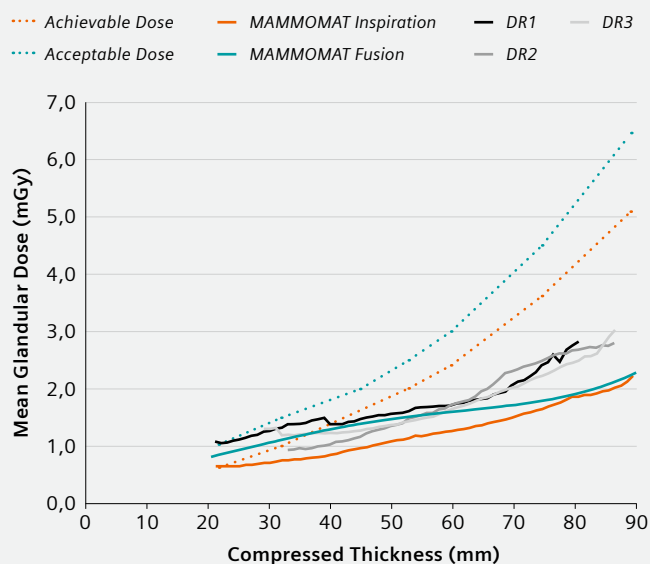

**Figure 1:** Average mean glandular dose as a function of compressed breast thickness. DR1 to DR3 are FFDM systems from different manufacturers.

## New generation CsI detector

The new generation Cesium-Iodide (CsI) high-resolution detector improves both spatial resolution and dose efficiency. This new detector technology is also remarkably robust, making MAMMOMAT Fusion particularly suitable for challenging situations, such as truck-installed mobile screening units.

The X-ray-sensitive area of each pixel of the new detector has been significantly increased with an innovative multi-layer design. This enables a smaller, 83 µm pixel size and higher detective quantum efficiency (DQE) compared to previous technologies. In addition, the scintillator thickness is matched to the X-ray spectrum of MAMMOMAT Fusion. The combination of these developments results in better image quality, at the right dose.

## Intelligent dose reduction


The combination of multiple features enables the intelligent dose reduction without compromising image quality.

OpDose® determines the exposure parameters required for each breast individually depending on breast density and thickness. OpDose® works in every position and automatically excludes the pectoral muscle from its calculation.

OpComp® with SoftSpeed slows the compression as soon as the plate reaches the breast. The compression is increased only as long as the breast is soft and pliable, stopping automatically when the optimal level of compression is reached, improving image quality and patient comfort.

OpView image processing software allows the selection from a range of predefined image impressions and flavors to match individual clinical requirements and preferences.

\*MAMMOMAT Fusion is not commercially available in all countries. Due to regulatory reasons its future availability cannot be guaranteed. Please contact your local Siemens organization for further information.

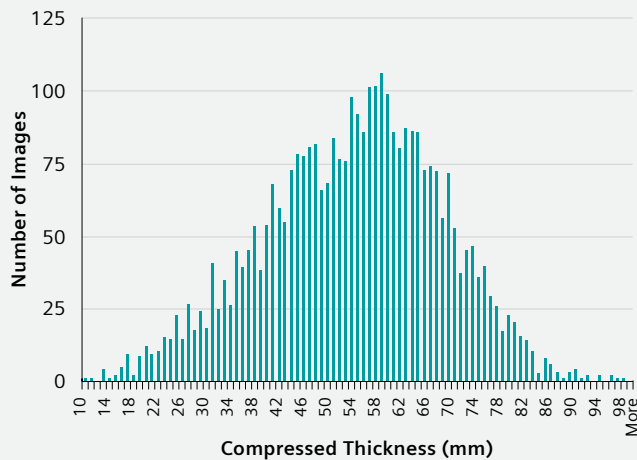

**Figure 2:** Histogram of 4000 images taken in a routine screening environment

## Conclusion

This study analyzed dose levels of clinical cases performed with Siemens MAMMOMAT Fusion, and compared them to dose levels of different FFDM systems obtained in a previous study<sup>1</sup>.

The mean glandular dose of MAMMOMAT Fusion is in the same range as other FFDM systems with a trend toward even lower dose with breast thickness above 50 mm.

MAMMOMAT Fusion's intelligent dose reduction features and new generation CsI detector deliver high image quality at the right dose.

## Dose comparison

Dose information was extracted from over 1000 cases (4000 images) performed with MAMMOMAT Fusion. The images were acquired and read in a routine screening environment under the appropriate image quality control conditions. A median dose value for each breast thickness was calculated allowing a dose – breast thickness curve to be plotted.

The result was then compared to the dose levels of a previous study 1, in order to obtain a comparison of the dose level of MAMMOMAT Fusion to other FFDM systems in the clinical routine. The resulting new dose – breast thickness curve for MAMMOMAT Fusion imposed on the curves from the previous study 1 can be seen in figure 1.

The curve for MAMMOMAT Fusion is in the same level range or below that of other FFDM systems. This is especially noticed at breast thicknesses higher than 50 mm, demonstrating the effectiveness of the intelligent dose reduction features and new generation Cesium-Iodide (CsI) high-resolution detector.

A histogram of the images taken (Figure 2) shows that the majority of examined breasts had more than 50 mm compressed thickness. Therefore, of the women examined in this sample, the majority have benefited from the lower dose with Siemens MAMMOMAT Fusion compared to other FFDM systems.

<sup>1</sup> Mean glandular dose survey of 2D mammograms acquired with the Siemens MAMMOMAT Inspiration system, White Paper, Hilde Bosmans, Koen Michielsens, Kim Lemmens, Jurgen Jacobs, 2011. Department of Radiology University, Hospital Leuven, Herestraat 49, B 3000 Leuven, Belgium.

Note: The above results do not allow a comparative evaluation of the intrinsic performance of the systems, as image quality factors or detector characteristics are not compared. Rather, it shows how the dose levels of these different DR modalities are preprogrammed on the modalities.

This document is not for distribution in the U.S. MAMMOMAT Fusion is not commercially available in all countries. Due to regulatory reasons its future availability cannot be guaranteed. Please contact your local Siemens Healthineers organisation for further information.

On account of certain regional limitations of sales rights and service availability, we cannot guarantee that all products/services/features included in this brochure are available through the Siemens Healthineers sales organization worldwide. Availability and packaging may vary by country and are subject to change without prior notice.

The information in this document contains general descriptions of the technical options available and may not always apply in individual cases. Siemens Healthineers reserves the right to modify the design and specifications contained herein without prior notice. Please contact your local Siemens Healthineers sales representative for the most current information.

In the interest of complying with legal requirements concerning the environmental compatibility of our products (protection of natural resources and waste conservation), we may recycle certain components where legally permissible. For recycled components we use the same extensive quality assurance measures as for factory-new components.

Any technical data contained in this document may vary within defined tolerances. Original images always lose a certain amount of detail when reproduced.

---

**Siemens Healthineers Headquarters**

Siemens Healthcare GmbH  
Henkestr. 127  
91052 Erlangen, Germany  
Phone: +49 9131 84-0  
[siemens-healthineers.com](http://siemens-healthineers.com)