



Release VA41A / 2024-10-02 / Revision 10242

syngo.share webview **syngo.share webview diagnostic**

DICOM 3.0 Conformance Statement

Overview

Content and transfer

Table 1: Supported Storage SOP Classes

SOP Class	SOP Class UID	DIMSE services		DICOM web services		Media services		Function			
		SCU	SCP	UA	OS	FSC	FSR	Create	Display	Process	Archive
Hardcopy Grayscale Image Storage	1.2.840.10008.5.1.1.29	Y	N	N	N	Y	N	N	Y	N	N
HardcopS Color Image Storage	1.2.840.10008.5.1.1.30	Y	N	N	N	Y	N	N	Y	N	N
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Y	N	N	N	Y	N	N	Y	N	N
Digital X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.1.1	Y	N	N	N	Y	N	N	Y	N	N
Digital X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.1.1.1	Y	N	N	N	Y	N	N	Y	N	N
Digital Mammography X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.1.2	Y	N	N	N	Y	N	N	Y	N	N
Digital Mammography X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.1.2.1	Y	N	N	N	Y	N	N	Y	N	N
Digital Intra-Oral X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.1.3	Y	N	N	N	Y	N	N	Y	N	N
Digital Intra-Oral X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.1.3.1	Y	N	N	N	Y	N	N	Y	N	N
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Y	N	N	N	Y	N	N	Y	Y	N
Enhanced CT Image Storage	1.2.840.10008.5.1.4.1.1.2.1	Y	N	N	N	Y	N	N	Y	N	N
Legacy Converted Enhanced CT Image Storage	1.2.840.10008.5.1.4.1.1.2.2	Y	N	N	N	Y	N	N	Y	N	N
Ultrasound Multi-frame Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.3	Y	N	N	N	Y	N	N	Y	N	N
Ultrasound Multi-frame Image Storage	1.2.840.10008.5.1.4.1.1.3.1	Y	N	N	N	Y	N	N	Y	N	N
MR Image Storage	1.2.840.10008.5.1.4.1.1.4	Y	N	N	N	Y	N	N	Y	Y	N
Enhanced MR Image Storage	1.2.840.10008.5.1.4.1.1.4.1	Y	N	N	N	Y	N	N	Y	N	N
Enhanced MR Color Image Storage	1.2.840.10008.5.1.4.1.1.4.3	Y	N	N	N	Y	N	N	Y	N	N
Legacy Converted Enhanced MR Image Storage	1.2.840.10008.5.1.4.1.1.4.4	Y	N	N	N	Y	N	N	Y	N	N
Nuclear Medicine Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.5	Y	N	N	N	Y	N	N	Y	N	N
Ultrasound Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.6	Y	N	N	N	Y	N	N	Y	N	N
Ultrasound Image Storage	1.2.840.10008.5.1.4.1.1.6.1	Y	N	N	N	Y	N	N	Y	N	N
Enhanced US Volume Storage	1.2.840.10008.5.1.4.1.1.6.2	Y	N	N	N	Y	N	N	Y	N	N
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	Y	N	N	N	Y	N	N	Y	N	N
Multi-frame Grayscale Byte Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.2	Y	N	N	N	Y	N	N	Y	N	N

Table 1 ▼

Table 1

SOP Class	SOP Class UID	DIMSE services		DICOM web services		Media services		Function			
		SCU	SCP	UA	OS	FSC	FSR	Create	Display	Process	Archive
Multi-frame Grayscale Word Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.3	Y	N	N	N	Y	N	N	Y	N	N
Multi-frame True Color Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.4	Y	N	N	N	Y	N	S	Y	N	N
12-lead ECG Waveform Storage	1.2.840.10008.5.1.4.1.1.9.1.1	Y	N	N	N	Y	N	N	Y	N	N
General ECG Waveform Storage	1.2.840.10008.5.1.4.1.1.9.1.2	Y	N	N	N	Y	N	N	Y	N	N
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1	Y	N	N	N	Y	N	N	Y	N	N
X-Ray Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.12.1	Y	N	N	N	Y	N	N	Y	N	N
Enhanced xA Image Storage	1.2.840.10008.5.1.4.1.1.12.1.1	Y	N	N	N	Y	N	N	Y	N	N
X-Ray Radiofluoroscopic Image Storage	1.2.840.10008.5.1.4.1.1.12.2	Y	N	N	N	Y	N	N	Y	N	N
Enhanced xRF Image Storage	1.2.840.10008.5.1.4.1.1.12.2.1	Y	N	N	N	Y	N	N	Y	N	N
X-Ray Angiographic Bi-Plane Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.12.3	Y	N	N	N	Y	N	N	Y	N	N
X-Ray 3D Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.13.1.1	Y	N	N	N	Y	N	N	Y	N	N
X-Ray 3D Craniofacial Image Storage	1.2.840.10008.5.1.4.1.1.13.1.2	Y	N	N	N	Y	N	N	Y	N	N
Breast Tomosynthesis Image Storage	1.2.840.10008.5.1.4.1.1.13.1.3	Y	N	N	N	Y	N	N	Y	N	N
Breast Projection X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.13.1.4	Y	N	N	N	Y	N	N	Y	N	N
Breast Projection X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.13.1.5	Y	N	N	N	Y	N	N	Y	N	N
Intravascular Optical Coherence Tomography Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.14.1	Y	N	N	N	Y	N	N	Y	N	N
Intravascular Optical Coherence Tomography Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.14.2	Y	N	N	N	Y	N	N	Y	N	N
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20	Y	N	N	N	Y	N	N	Y	N	N
VL Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.77.1	Y	N	N	N	Y	N	N	Y	N	N
VL Endoscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.1	Y	N	N	N	Y	N	N	Y	N	N
Video Endoscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.1.1	Y	N	N	N	Y	N	N	Y	N	N
VL Microscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.2	Y	N	N	N	Y	N	N	Y	N	N
Video Microscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.2.1	Y	N	N	N	Y	N	N	Y	N	N
VL Slide-Coordinates Microscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.3	Y	N	N	N	Y	N	N	Y	N	N
VL Photographic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.4	Y	N	N	N	Y	N	N	Y	N	N
Video Photographic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.4.1	Y	N	N	N	Y	N	N	Y	N	N
Ophthalmic Photography 8 Bit Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.1	Y	N	N	N	Y	N	N	Y	N	N
Ophthalmic Photography 16 Bit Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.2	Y	N	N	N	Y	N	N	Y	N	N
Ophthalmic Tomography Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.4	Y	N	N	N	Y	N	N	Y	N	N

Table 1

^ Table 1

SOP Class	SOP Class UID	DIMSE services		DICOM web services		Media services		Function			
		SCU	SCP	UA	OS	FSC	FSR	Create	Display	Process	Archive
Wide Field Ophthalmic Photography Stereographic Projection Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.5	Y	N	N	N	Y	N	N	Y	N	N
Wide Field Ophthalmic Photography 3D Coordinates Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.6	Y	N	N	N	Y	N	N	Y	N	N
VL Whole Slide Microscopy Image Storage	1.2.840.10008.5.1.4.1.1.77.1.6	Y	N	N	N	Y	N	N	Y	N	N
VL Multi-frame Image Storage - Trial (Retired)	1.2.840.10008.5.1.4.1.1.77.2	Y	N	N	N	Y	N	N	Y	N	N
Ophthalmic Thickness Map Storage	1.2.840.10008.5.1.4.1.1.81.1	Y	N	N	N	Y	N	N	Y	N	N
Corneal Topography Map Storage	1.2.840.10008.5.1.4.1.1.82.1	Y	N	N	N	Y	N	N	Y	N	N
Basic Text SR	1.2.840.10008.5.1.4.1.1.88.11	Y	N	N	N	Y	N	N	Y	N	N
Key Object Selection Document Storage	1.2.840.10008.5.1.4.1.1.88.59	Y	N	N	N	Y	N	N	Y	N	N
X-Ray Radiation Dose SR Storage	1.2.840.10008.5.1.4.1.1.88.67	Y	N	N	N	Y	N	N	Y	N	N
Radiopharmaceutical Radiation Dose SR Storage	1.2.840.10008.5.1.4.1.1.88.68	Y	N	N	N	Y	N	N	Y	N	N
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1	Y	N	N	N	Y	N	N	Y	N	N
Positron Emission Tomography Image Storage	1.2.840.10008.5.1.4.1.1.128	Y	N	N	N	Y	N	N	Y	Y	N
Legacy Converted Enhanced PET Image Storage	1.2.840.10008.5.1.4.1.1.128.1	Y	N	N	N	Y	N	N	Y	N	N
Enhanced PET Image Storage	1.2.840.10008.5.1.4.1.1.130	Y	N	N	N	Y	N	N	Y	N	N
RT Image Storage	1.2.840.10008.5.1.4.1.1.481.1	Y	N	N	N	Y	N	N	Y	N	N
RT Dose Storage	1.2.840.10008.5.1.4.1.1.481.2	Y	N	N	N	Y	N	N	Y	N	N
DICOS CT Image Storage	1.2.840.10008.5.1.4.1.1.501.1	Y	N	N	N	Y	N	N	Y	N	N

Table 2: Supported Storage Transfer Syntaxes

Transfer Syntax set	Transfer Syntax name	Transfer Syntax UID
Lossless Compressed Transfer Syntax Set (LL)	Deflated Explicit VR Little Endian	1.2.840.10008.1.2.1.99
	JPEG Lossless, Non-Hierarchical (Process 14)	1.2.840.10008.1.2.4.57
	JPEG Lossless, Non-Hierarchical, First-Order Prediction (Process 14 [Selection Value 1]): Default Transfer Syntax for Lossless JPEG Image Compression	1.2.840.10008.1.2.4.70
	JPEG-LS Lossless Image Compression	1.2.840.10008.1.2.4.80
	JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90
	RLE Lossless	1.2.840.10008.1.2.5
Lossy Compressed Transfer Syntax Set (L)	JPEG Baseline (Process 1): Default Transfer Syntax for Lossy JPEG 8-bit Image Compression	1.2.840.10008.1.2.4.50
	JPEG Extended (Processes 2 & 4): Default Transfer Syntax for Lossy JPEG 12-bit Image Compression (Process 4 only)	1.2.840.10008.1.2.4.51
	JPEG-LS Lossy (Near-Lossless) Image Compression	1.2.840.10008.1.2.4.81
	JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91
Non-Image Transfer Syntax Set (NI)	Implicit VR Little Endian: Default Transfer Syntax for DICOM	1.2.840.10008.1.2
	Explicit VR Little Endian	1.2.840.10008.1.2.1
	Deflated Explicit VR Little Endian	1.2.840.10008.1.2.1.99
	Explicit VR Big Endian	1.2.840.10008.1.2.2
Uncompressed Transfer Syntax Set (U)	Implicit VR Little Endian: Default Transfer Syntax for DICOM	1.2.840.10008.1.2
	Explicit VR Little Endian	1.2.840.10008.1.2.1
	Explicit VR Big Endian	1.2.840.10008.1.2.2
Video Transfer Syntax Set (V)	MPEG2 Main Profile / Main Level	1.2.840.10008.1.2.4.100
	MPEG2 Main Profile / High Level	1.2.840.10008.1.2.4.101
	MPEG-4 AVC/H.264 High Profile / Level 4.1	1.2.840.10008.1.2.4.102
	MPEG-4 AVC/H.264 BD-compatible High Profile / Level 4.1	1.2.840.10008.1.2.4.103
	MPEG-4 AVC/H.264 High Profile / Level 4.2 For 2D Video	1.2.840.10008.1.2.4.104
	MPEG-4 AVC/H.264 High Profile / Level 4.2 For 3D Video	1.2.840.10008.1.2.4.105
	MPEG-4 AVC/H.264 Stereo High Profile / Level 4.2	1.2.840.10008.1.2.4.106
	HEVC/H.265 Main Profile / Level 5.1	1.2.840.10008.1.2.4.107
	HEVC/H.265 Main 10 Profile / Level 5.1	1.2.840.10008.1.2.4.108

DIMSE services

Query/Retrieve

Table 3: Supported Query/Retrieve SOP Classes

SOP Class	SOP Class UID	Transfer Syntax	Transfer Syntax UID	SCU	SCP
Study Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	Implicit VR Little Endian	1.2.840.10008.1.2	Y	N
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Study Root Query/Retrieve Information Model – MOVE	1.2.840.10008.5.1.4.1.2.2.2	Implicit VR Little Endian	1.2.840.10008.1.2	Y	N
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

Contents

Overview	2
Content and transfer	2
DIMSE services	6
Query/Retrieve	6
1 Introduction	10
1.1 Remarks	10
2 Implementation model	11
2.1 Application data flow diagram	11
2.2 Functional definition of AEs	11
2.3 Sequencing of real-world activities	12
3 Application Entity specifications	13
3.1 Supported SOP Classes and Transfer Syntaxes	13
3.2 Association establishment policies	13
3.2.1 General	13
3.2.2 Number of associations	13
3.2.3 Asynchronous nature	13
3.2.4 Implementation identifying information	13
3.2.5 Association initiation policy by real-world activity	13
4 DICOM media AE specification	14
4.1 Implementation model	14
4.1.1 Application data flow diagram	14
4.1.2 Functional definitions of the AE	14
4.1.3 Sequencing of real-world activities	14
4.1.4 File meta information	14
4.2 Application Entity specification	14
4.2.1 Real-world activities	14
4.3 Augmented and private application profiles	15
5 Communication profiles	16
5.1 Supported communication stacks	16
6 Configuration	17
7 Support of extended character sets	18
7.1 Supported character sets	18
7.2 Usage of Specific Character Set in C-FIND and C-MOVE requests	18
A DICOM attribute list for Query/Retrieve service classes	19
A.1 Used C-FIND attributes in requests issued by <i>syngo.share</i> webview	19
B Limitations	21
C webrenderserver	24
C.1 Implementation model	24
C.1.1 <i>webrenderserver</i> data flow diagram	24
C.1.2 Functional definition of AEs	24
C.1.3 Sequencing of real-world activities	24

C.2	Application Entity specifications	24
C.2.1	WADO-RS general notes	24
C.2.2	WADO-RS rendered resources	25
C.2.3	WADO-RS thumbnail resources	26
C.2.4	WADO-RS supported Transfer Syntaxes	27
C.2.5	WADO-RS connection policies	27

List of figures

1	<i>syngo.share</i> webview FSC data flow diagram	11
2	<i>syngo.share</i> webview AE application data flow diagram	11
3	<i>webrenderserver</i> data flow diagram	24

List of tables

1	Supported Storage SOP Classes	2
2	Supported Storage Transfer Syntaxes	5
3	Supported Query/Retrieve SOP Classes	6
4	<i>syngo.share</i> webview AE-related application profiles, real-world activity roles, and roles for interchanging	14
5	Supported <i>syngo.share</i> webview Export Transfer Syntaxes	15
6	Supported character sets	18
7	Used C-FIND attributes in requests issued by <i>syngo.share</i> webview AE (C-FIND SCU) on the study level	19
8	Used C-FIND attributes in requests issued by <i>syngo.share</i> webview AE (C-FIND SCU) on the series level	19
9	Used C-FIND attributes in requests issued by <i>syngo.share</i> webview AE (C-FIND SCU) on the image level	20
10	SOP Classes and corresponding UIDs for which Grayscale Softcopy Presentation States are not applied	22
11	Grayscale Softcopy Presentation State module coverage and limitations	23
12	Rendered resources	25
13	Options and restrictions on rendered resources	25
14	Thumbnail resources	26
15	Options and restrictions on thumbnail resources	26
16	WADO-RS supported Transfer Syntaxes	27
17	WADO-RS HTTP status codes	28

1 Introduction

This document is a DICOM 3.0 Conformance Statement that describes the DICOM capabilities of *syngo.share webview*.

syngo.share webview is a versatile multi-modality display system for DICOM images. It is able to retrieve and display DICOM images from *syngo.share core* or, via Query/Retrieve, from third-party DICOM archives. Additionally, *syngo.share webview* supports the export of DICOM images, series and studies.

1.1 Remarks

This Conformance Statement is intended to aid in the validation of the integration of *syngo.share webview* within a DICOM environment. This statement is not meant to replace the validation with other DICOM equipment to ensure the intended, proper exchange of information. Thus, it is still important to ensure the proper interoperability of the intended DICOM integration.

The user must be aware of the following issues:


- The comparison of different Conformance Statements should be the first step towards an assessment of the interoperability within a DICOM environment.
- Testing procedures should be defined to validate the desired level of connectivity.

2 Implementation model

2.1 Application data flow diagram


syngo.share webview provides a user interface for exporting DICOM files – the FSC (File Set Creator) – to portable media (such as USB sticks), network directories or the local file system. These functions are integral components of *syngo.share* webview. Within the context of this Conformance Statement, the exporting functionality is referred to as the *syngo.share* webview FSC AE.

Figure 1: *syngo.share* webview FSC data flow diagram

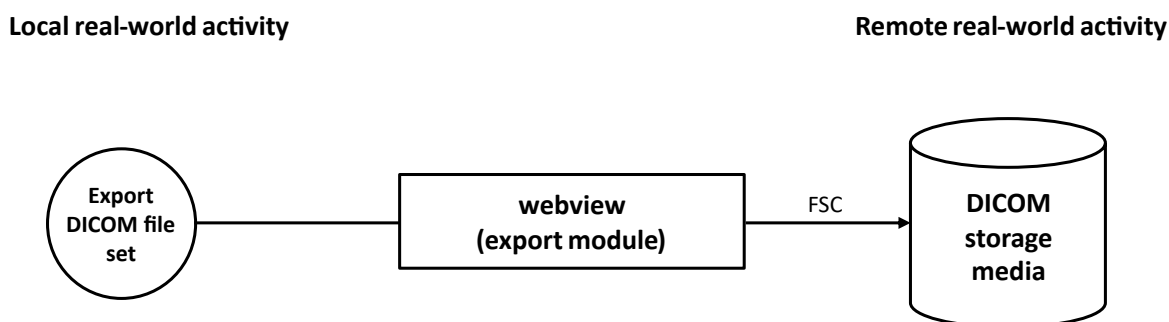

Additionally, *syngo.share* webview provides the ability to query third-party PACS and retrieve data from them. *syngo.share* webview also provides report-reading capabilities for Structured Reports and the possibility to apply Grayscale Softcopy Presentation States (GSPS) as defined in DICOM PS3.3.

Figure 2: *syngo.share* webview AE application data flow diagram

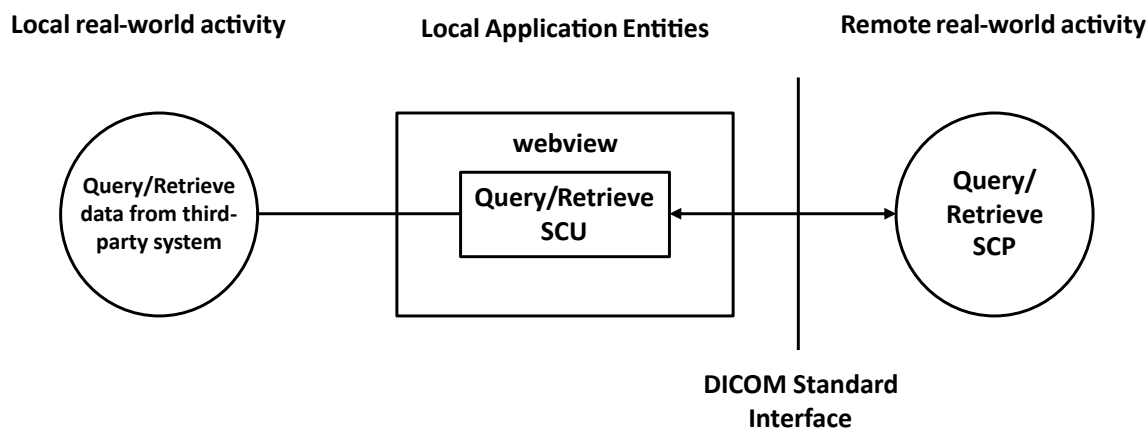

2.2 Functional definition of AEs

The *syngo.share* webview FSC AE is able to write user-selected DICOM files like DICOMDIRs or image objects which are compliant with DICOM PS3.10. These files can be located either on the local file system or on DICOM PS3.12-compliant media.

The *syngo.share* webview AE is able to query third-party PACS using the C-FIND service and retrieve data using the C-MOVE service in accordance with DICOM PS3.7. The third-party PACS can be queried

using the Study Root Query/Retrieve Information Model as defined in DICOM PS3.4, C.3.2. Baseline behavior is supported for C-FIND SCU and C-MOVE SCU.

The *syngo.share* webview AE is able to render DICOM instances received via Query/Retrieve. Section 3.1 describes the supported SOP Classes for *syngo.share* webview.

2.3 Sequencing of real-world activities

- The user can request the services of the *syngo.share* webview AE at any time through the user interface.

3 Application Entity specifications

The *syngo.share* webview FSC provides functionalities for the sole purpose of handling DICOM media. The DICOM Storage SOP Classes listed in Table 1 are supported for viewing. For SOP Classes not listed in Table 1, accessing the DICOM header is possible within the application.

3.1 Supported SOP Classes and Transfer Syntaxes

The supported DICOM SOP Classes are described in Table 1 and Table 3. Table 2 lists the supported storage Transfer Syntaxes.

Note that SOP Classes for the purpose of image data storage require the DICOM Image Pixel module to exist in order to display image data. Also, a value must be provided in the DICOM element Pixel Data (7FE0, 0010), and the underlying Transfer Syntax must be supported. For limitations please see Appendix B.

3.2 Association establishment policies

3.2.1 General

The *syngo.share* webview AE supports TCP/IP. When a user requests a C-FIND or C-MOVE operation, it attempts to establish an association with a remote AE. The host, port and remote AE title are defined within the server configuration. The maximum accepted PDU size is 16378.

3.2.2 Number of associations

The *syngo.share* webview AE supports multiple simultaneous associations for C-FIND, C-MOVE and C-STORE operations. The *syngo.share* webview AE is configured for ten simultaneous connections for C-STORE by default. This value can be changed but requires a restart.

3.2.3 Asynchronous nature

The *syngo.share* webview AE only allows a single outstanding operation on each association, meaning that it does not perform asynchronous negotiation.

3.2.4 Implementation identifying information

- **Implementation Class UID**
1.2.40.0.13.1.3
- **Implementation Version Name**
dcm4che-5.31.2

3.2.5 Association initiation policy by real-world activity

The *syngo.share* webview AE initiates an association with a remote AE for C-FIND and C-MOVE requests. The proposed Transfer Syntaxes for association establishment for C-FIND and C-MOVE are defined in Table 3.

4 DICOM media AE specification

This chapter describes the DICOM media functionalities of the *syngo.share* webview AE.

4.1 Implementation model

4.1.1 Application data flow diagram

See Section 2.1 (*syngo.share* webview data flow diagram).

4.1.2 Functional definitions of the AE

The *syngo.share* webview AE implements standard DICOM-conformant service classes for the creation of DICOM file sets (according to DICOM PS3.10).

4.1.3 Sequencing of real-world activities

The DICOM media functionalities of the *syngo.share* webview AE can be used at any time through its user interface.

4.1.4 File meta information

- **Implementation Class UID**
1.2.40.0.13.1.3
- **Implementation Version Name**
dcm4che-5.31.2

4.2 Application Entity specification

See Section 3 for supported SOP Classes for the export of media.

Table 4: *syngo.share* webview AE-related application profiles, real-world activity roles, and roles for interchanging

Supported application profiles	Real-world activity	Role	SC option
<i>syngo.share</i> webview FSC	Export DICOM file set	FSC	Interchange

4.2.1 Real-world activities

4.2.1.1 Creating DICOMDIRs

The *syngo.share* webview AE supports the creation of DICOMDIR with all mandatory keys, as defined in DICOM PS3.10.

Exporting media

The *syngo.share* webview AE is able to organize DICOM images, series and studies into a file set which is then saved to the local file system or to DICOM PS3.12-compliant media. The *syngo.share* webview AE, here acting as an FSC, uses the following Transfer Syntaxes:

Table 5: Supported *syngo.share* webview Export Transfer Syntaxes

Transfer Syntax name	Transfer Syntax UID
Explicit VR Little Endian	1.2.840.10008.1.2.1

The *syngo.share* webview AE also supports the usage of the original Transfer Syntax if the DICOM instance cannot be converted to the Explicit VR Little Endian Syntax. In this case, the created DICOMDIR is not compliant with the DICOM Standard. This option can be deactivated in *syngo.share* webview.

4.3 Augmented and private application profiles

Not used.

5 Communication profiles

5.1 Supported communication stacks

`syngo.share webview` provides plain TCP (see DICOM PS3.8, 9) communication. It uses the `dcm4che` library and toolkit for its communication.

6 Configuration

`syngo.share` webview is configured in `syngo.share` Configuration and static configuration files. See the product-specific technical documentation for further information.

7 Support of extended character sets

7.1 Supported character sets

Table 6 contains the character sets which are supported with and without code extension techniques.

Table 6: Supported character sets

MIME name	Without code extensions	With code extensions
US-ASCII		ISO 2022 IR 6
ISO-8859-1	ISO_IR 100	ISO 2022 IR 100
ISO-8859-2	ISO_IR 101	ISO 2022 IR 101
ISO-8859-3	ISO_IR 109	ISO 2022 IR 109
ISO-8859-4	ISO_IR 110	ISO 2022 IR 110
ISO-8859-5	ISO_IR 144	ISO 2022 IR 144
ISO-8859-6	ISO_IR 127	ISO 2022 IR 127
ISO-8859-7	ISO_IR 126	ISO 2022 IR 126
ISO-8859-8	ISO_IR 138	ISO 2022 IR 138
ISO-8859-9	ISO_IR 148	ISO 2022 IR 148
JIS_X0201	ISO_IR 13	ISO 2022 IR 13
TIS-620	ISO_IR 166	ISO 2022 IR 166
JIS_X0208-1990	-	ISO 2022 IR 87
JIS_X0212-1990	-	ISO 2022 IR 159
KS_X_1001-1997	-	ISO 2022 IR 149
GB2312	-	ISO 2022 IR 58
UTF-8	ISO_IR 192	-
GB18030	GB18030	-
GBK	GBK	-

7.2 Usage of Specific Character Set in C-FIND and C-MOVE requests

syngo.share webview AE as C-FIND and C-MOVE SCU uses Specific Character Set (0008,0005) with a value of ISO_IR 192 for requests.

A DICOM attribute list for Query/Retrieve service classes

The following tables list the DICOM keys used for matching and/or return on the study (Table 7), series (Table 8) and image (Table 9) level in C-FIND requests issued by *syngo.share* webview AE as C-FIND SCU.

The Query Retrieve Level (0008,0052) and Specific Character Set (0008,0005) attributes are set in every C-FIND request and are therefore not mentioned in the tables below. The value of Query Retrieve Level is set to the respective level (study, series, image). The value of Specific Character Set is set to ISO_IR 192. In addition, the Unique Key Attribute for each level above the Query/Retrieve level is sent as well with a single value (baseline behavior).

A.1 Used C-FIND attributes in requests issued by *syngo.share* webview

Table 7: Used C-FIND attributes in requests issued by *syngo.share* webview AE (C-FIND SCU) on the study level

Attribute name	Tag	Usage
StudyDate	(0008,0020)	Matching and Return
StudyTime	(0008,0030)	Matching and Return
AccessionNumber	(0008,0050)	Matching and Return
ModalitiesInStudy	(0008,0061)	Matching and Return
ReferringPhysicianName	(0008,0090)	Matching and Return
StudyDescription	(0008,1030)	Matching and Return
PatientName	(0010,0010)	Matching and Return
PatientID	(0010,0020)	Matching and Return
PatientBirthDate	(0010,0030)	Matching and Return
PatientBirthTime	(0010,0032)	Return
PatientSex	(0010,0040)	Matching and Return
PatientBirthName	(0010,1005)	Matching and Return
PatientMotherBirthName	(0010,1060)	Matching and Return
StudyInstanceUID	(0020,000D)	Matching and Return
StudyID	(0020,0010)	Matching and Return
NumberOfStudyRelatedSeries	(0020,1206)	Return
NumberOfStudyRelatedInstances	(0020,1208)	Return
AdmissionID	(0038,0010)	Return

Table 8: Used C-FIND attributes in requests issued by *syngo.share* webview AE (C-FIND SCU) on the series level

Attribute name	Tag	Usage
SeriesDate	(0008,0021)	Return
SeriesTime	(0008,0031)	Return
Modality	(0008,0060)	Return
Manufacturer	(0008,0070)	Return

Table 8 ▼

^ Table 8

Attribute name	Tag	Usage
InstitutionName	(0008,0080)	Matching and Return
StationName	(0008,1010)	Matching and Return
SeriesDescription	(0008,103E)	Return
InstitutionalDepartmentName	(0008,1040)	Matching and Return
PerformingPhysicianName	(0008,1050)	Return
OperatorsName	(0008,1070)	Return
ManufacturerModelName	(0008,1090)	Return
BodyPartExamined	(0018,0015)	Matching and Return
SeriesInstanceUID	(0020,000E)	Matching and Return
SeriesNumber	(0020,0011)	Matching and Return
Laterality	(0020,0060)	Matching and Return
NumberOfSeriesRelatedInstances	(0020,1209)	Return
PerformedProcedureStepStartDate	(0040,0244)	Matching and Return
PerformedProcedureStepStartTime	(0040,0245)	Return
RequestAttributesSequence/ ScheduledProcedureStepID	(0040,0275)/(0040,0009)	Matching and Return
RequestAttributesSequence/ RequestedProcedureID	(0040,0275)/(0040,1001)	Matching and Return

Table 9: Used C-FIND attributes in requests issued by syngo.share webview AE (C-FIND SCU) on the image level

Attribute name	Tag	Usage
ImageType	(0008,0008)	Return
SOPClassUID	(0008,0016)	Return
SliceThickness	(0018,0050)	Return
AcquisitionNumber	(0020,0012)	Return
InstanceNumber	(0020,0013)	Return
ImagePositionPatient	(0020,0032)	Return
ImageOrientationPatient	(0020,0037)	Return
SliceLocation	(0020,1041)	Return
SamplesPerPixel	(0028,0002)	Return
PhotometricInterpretation	(0028,0004)	Return
NumberOfFrames	(0028,0008)	Return
Rows	(0028,0010)	Return
Columns	(0028,0011)	Return
BitsAllocated	(0028,0100)	Return
BitsStored	(0028,0101)	Return
ConceptNameCodeSequence/ CodeValue	(0040,A043)/(0008,0100)	Matching and Return
ConceptNameCodeSequence/ CodingSchemeDesignator	(0040,A043)/(0008,0102)	Matching and Return
ConceptNameCodeSequence/ CodingSchemeVersion	(0040,A043)/(0008,0103)	Return
ConceptNameCodeSequence/ CodeMeaning	(0040,A043)/(0008,0104)	Return
ReferencedRequestSequence/ AccessionNumber	(0040,A370)/(0008,0050)	Return
ReferencedRequestSequence/ StudyInstanceUID	(0040,A370)/(0020,000D)	Return
ReferencedRequestSequence/ RequestedProcedureID	(0040,A370)/(0040,1001)	Return
ContentLabel	(0070,0080)	Return
ContentDescription	(0070,0081)	Return
PresentationCreationDate	(0070,0082)	Return
PresentationCreationTime	(0070,0083)	Return
ContentCreatorName	(0070,0084)	Return

B Limitations

DICOM instances

- Timezone Offset From UTC (0008,0201) is not applied on DA, TM, and DT attributes.
- Instances for which the value length of a tag exceeds 2 GiB (2^{31} - 9 bytes) are not supported.

DICOM images

- Images with 1 bit allocated for each pixel sample (binary images) are not supported.
- Images with more than 16 bits per sample (Bits Stored (0028,0101)) are not supported.
- Multiple window center and width presets are not supported.
- VOI LUTs present in the Functional Group of multiframe images are not considered.
- VOI LUT functions are not supported.
- Overlays stored in the unused bit planes of Pixel Data (7FE0,0010) are not supported, since this method was retired from the DICOM Standard in 2004.
- Display Shutter and other graphical modules defined in image IODs are not applied.
- Length calculations on ultrasound images are not supported.
- ECG visualization is available only for color and monochrome images of modality XA. Only the first repeating group (5000 group) of the Curve Module is considered if Type of Data (5000,0020) is ECG. In addition, Curve Dimensions (5000,0005) must be 2, and Curve Data Descriptor (5000,0110) must be present with values 01. Value representation is limited to unsigned short data: Data Value Representation (5000,0103) must be 0 and VR of Curve Data (5000,3000) must be OW.

DICOM whole slide microscopy images

- Images are interpreted as common multiframe images if DICOM header- and frame-based access is not available.
- Monochrome images are currently not supported.
- Images for which the Dimension Organization Type (0020,9311) is present with value must either follow TILED_FULL or TILED_SPARSE dimension organization.
- All images with an Image Flavor (third value of Image Type (0008,0008)) of VOLUME must use the same dimension organization.
- Thumbnail and overview layer are not available if the edge length of the topmost image of the image pyramid exceeds 5,000 pixels.
- Overview layer is not available if the edge length of the topmost image of the image pyramid falls below 800 pixels.
- Multiple items in the Dimension Organization Sequence (0020,9221) and Dimension Index Sequence (0020A,9222) are not supported.
- Images with multiple optical paths are not supported.
- Images with Dimension Index Pointer (003A,0210) not targeting the Plane Position Slide Sequence (0048,021A) are not supported.

- All images of the underlying series must share the same Frame of Reference UID (0020,0052), Container Identifier (0040,0512), Issuer of the Container Identifier Sequence (0040,0513), Specimen UID (0040,0554), and Image Orientation (Slide) (0048,0102).
- Images with an Image Flavor (third value of Image Type (0008,0008)) of LABEL or OVERVIEW consisting of multiple frames are not supported.
- The number of Z-layers and the spacing between slices must be uniform across all zoom layers.
- Images with a layer width or height exceeding $2^{31} - 1$ pixels are not supported.
- Images with a number of Z-layers exceeding $2^{31} - 1$ are not supported.
- Images with layers exceeding $2^{31} - 9$ tiles (frames) are not supported.
- Images with Dimension Organization Type (0020,9311) of TILED_FULL and with an Image Flavor (third value of Image Type (0008,0008)) of VOLUME are not supported if they are part of a concatenation or if the number of frames per layer or in total exceeds $2^{31} - 1$.

DICOM ECGs

- Only ECG SOP Classes listed in Table 1 are currently supported.
- Channel Sensitivity (003A,0210) is required to be present since samples shall represent defined (not arbitrary) units.
- An ECG is not supported if Channel Time Skew (003A,0214), Channel Sample Skew (003A,0215), or Channel Offset (003A,0218) is present with a non-zero value.
- Channel layouts are only possible if exactly 12 channels are present in the selected Waveform Multiplex Group.
- Display properties (e.g. presentation groups, display colors, display scale) of the Waveform Module are not supported.

DICOM videos

- Videos with an encapsulated video data stream segmented into more than one fragment are not supported.

DICOM Grayscale Softcopy Presentation States

- GSPS functionality is limited to viewing, as no GSPS are created.
- GSPS are not applied on instances of the SOP Classes listed in Table 10 due to their limited support.
- GSPS module limitations are listed in Table 11.
- Only the Presentation State Storage SOP Classes listed above (see Table 1) are currently supported.

Table 10: SOP Classes and corresponding UIDs for which Grayscale Softcopy Presentation States are not applied

SOP Class name	SOP Class UID
Segmentation Storage	1.2.840.10008.5.1.4.1.1.66.4
Surface Segmentation Storage	1.2.840.10008.5.1.4.1.1.66.5
VL Whole Slide Microscopy Image Storage	1.2.840.10008.5.1.4.1.1.77.1.6

Table 11: Grayscale Softcopy Presentation State module coverage and limitations

Module	Implemented	Limitations
Presentation State Identification	Partially	No further attributes beside Presentation Creation Date (0070,0082), Presentation Creation Time (0070,0083), and Content Label (0070,0080) with Content Creator's Name (0070,0084) from the Content Identification Macro are processed.
Presentation State Relationship	Partially	The Referenced Segment Number (0062,000B) from the Image SOP Instance Reference Macro is not considered.
Presentation State Shutter	No	-
Presentation State Mask	No	-
Mask	No	-
Display Shutter	No	-
Bitmap Display Shutter	No	-
Overlay Plane	No	-
Overlay Activation	No	-
Displayed Area	Partially	The Referenced Segment Number (0062,000B) from the Image SOP Instance Reference Macro is not considered. The enumerated value VOLUME of Pixel Origin Interpretation (0048,0301) is not supported. The enumerated value TRUE SIZE of Presentation Size Mode (0048,0301) is interpreted as a magnification factor of 1: 1, as the client currently does not provide the physical pixel size of the display.
Graphic Annotation	Partially	The Referenced Segment Number (0062,000B) from the Image SOP Instance Reference Macro is not considered. The enumerated value MATRIX for all annotation units of this module is not supported. The Text Style Sequence Macro, Compound Graphic Instance ID (0070,0226), and Graphic Group ID (0070,0295) from the Text Object Sequence (0070,0008) are not supported. The Line Style Sequence Macro, Fill Style Sequence Macro, Compound Graphic Instance ID (0070,0226), and Graphic Group ID (0070,0295) from the Graphic Object Sequence are not supported. The Compound Graphic Sequence is not considered.
Spatial Transformation	Yes	-
Graphic Layer	Partially	The Graphic Layer Description (0070,0068) is not considered.
Graphic Group	No	-
Modality LUT	No	-
Softcopy VOI LUT	No	-
Softcopy Presentation LUT	No	-
SOP Common	Partially	No further information is read from this module for GSPS instances apart from common information required for every supported DICOM instance (e.g. Specific Character Set (0008,0005)).

DICOM Structured Reports

- Support is currently limited to the rendering of Basic Text Structured Reports, X-Ray Radiation Dose Structured Reports, Radiopharmaceutical Radiation Dose Structured Reports, and Key Object Selection Documents.
- Referenced images need to be in the same study in order to be resolved automatically.
- Referenced images do not consider a referenced Presentation State or Segment Number.
- Referenced multiframe images not referencing all frames only consider the first referenced frame.


C webrenderserver

C.1 Implementation model

webrenderserver provides a WADO-RS interface for the retrieval of rendered and thumbnail resources. It receives and responds to WADO-RS requests from a remote AE targeting on DICOM images stored in *syngo.share* core.

C.1.1 *webrenderserver* data flow diagram

Figure 3: *webrenderserver* data flow diagram

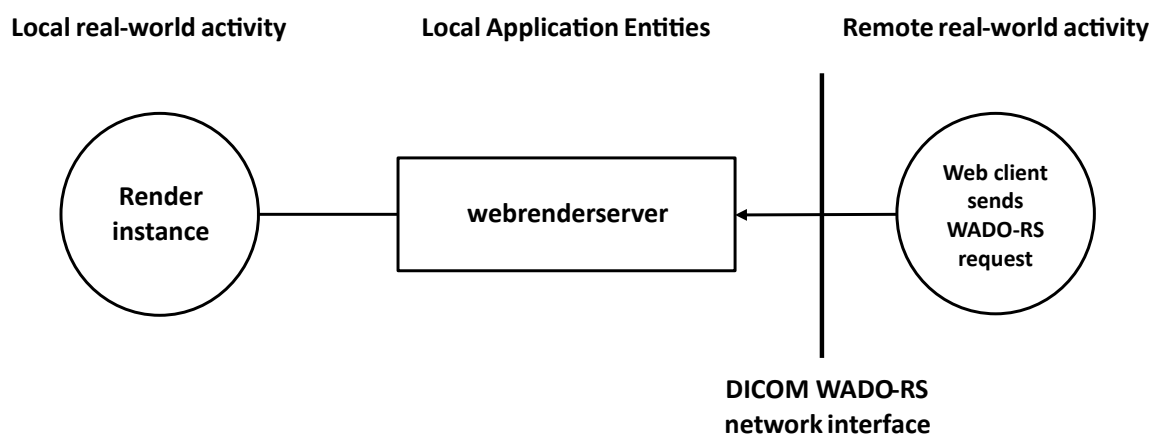

C.1.2 Functional definition of AEs

The *webrenderserver* AE waits for incoming WADO-RS requests from remote AEs (e.g. from a web browser) for rendered and thumbnail resources according to DICOM PS3.18. A remote AE can request rendered and thumbnail representations of DICOM images stored in *syngo.share* core on study, series, instance, and frame level.

C.1.3 Sequencing of real-world activities

- The user or the application performing the WADO-RS request can also request the services of the *webrenderserver* AE at any time.

C.2 Application Entity specifications

The *webrenderserver* AE provides a WADO-RS interface for the retrieval of rendered and thumbnail resources, as defined in and implemented according to DICOM PS3.18.

C.2.1 WADO-RS general notes

A remote AE can request rendered and thumbnail representations of DICOM images (color and grayscale) stored in *syngo.share* core on study, series, instance, and frame level. WADO-RS requests are syntactically checked and are required to contain valid authorization information by means of URL

signatures (HMAC) before getting processed. DICOM data for the requested resource are queried and are loaded from *syngo.share* core, rendered according to the specified options in the request, and sent to the remote AE as a single- or multipart HTTP response. The HTTP status code of the WADO-RS response informs the remote AE about the outcome of processing the WADO-RS request.



Note

The WADO-RS interface supports the retrieve transaction on rendered resources as well as on thumbnail resources. Further target resources (e.g. DICOM resources) are not supported. Furthermore, the retrieve capabilities transaction is not implemented due to security considerations (WADL disclosure).

C.2.2 WADO-RS rendered resources

The supported rendered resources along with their associated URI templates are listed in Table 12.

Table 12: Rendered resources

Resource	URI template
Rendered Study	/studies/{study}/rendered
Rendered Series	/studies/{study}/series/{series}/rendered
Rendered Instance	/studies/{study}/series/{series}/instances/{instance}/rendered
Rendered Frames	/studies/{study}/series/{series}/instances/{instance}/frames/{frames}/rendered

The options and restrictions valid for all rendered resources are listed in Table 13.

Table 13: Options and restrictions on rendered resources

Options	Restrictions
Accept (Request Header Field)	Restricted to image/jpeg, image/gif, image/png, or wildcard representations covering the three rendered media types.
Accept-Charset (Request Header Field)	Not applied as the response payload consists of bitmap images only.
Accept (Query Parameter)	Restricted to image/jpeg, image/gif, or image/png.
Charset (Query Parameter)	Not applied as the response payload consists of bitmap images only.
Annotation (Query Parameter)	Restricted to patient and/or technique. Other keywords are ignored and listed in a warning header in the response. Localization of the burned-in annotations (e.g. regarding date format) is based on the Accept-Language header field.
Quality (Query Parameter)	Restricted to JPEG quality if image/jpeg is requested.
Viewport (Query Parameter)	Restricted to vw (viewport width) and vh (viewport height). Further parameter values are ignored. If viewport width and height are specified, the rendered images will be of the specified dimensions with black background color if padding is necessary.
Transfer Syntaxes Supported	Only DICOM instances of a supported Transfer Syntax are considered for rendering. For the list of supported Transfer Syntaxes, see Table 16.
SOP Class Restrictions	Only DICOM instances of a supported SOP Class are considered for rendering. Rendering is restricted to SOP Classes for the purpose of image data storage - e.g. the DICOM Image Pixel Module exists along with the DICOM Element PixelData (7FE0, 0010).
Image Type Restrictions	Binary images are not supported.
Size Restriction	There are no size limits imposed by the WADO-RS AE. However, depending on system resource and file format restrictions there might be limits to the size of each image.

Table 13 ▼

^ Table 13

Options	Restrictions
Multi-Frame Image Handling	Each frame of a multi-frame image is rendered separately. Hence, each frame is returned as a part in case of a multi-part response. This is also the case if image/gif is the requested rendered media type.
Frames	The given frame numbers for rendered frames are expected to be in ascending order. If this is not the case, automatic sorting applies.

C.2.3 WADO-RS thumbnail resources

The supported thumbnail resources along with their associated URI templates are listed in Table 14.

Table 14: Thumbnail resources

Resource	URI template
Study Thumbnail	/studies/{study}/thumbnail
Series Thumbnail	/studies/{study}/series/{series}/thumbnail
Instance Thumbnail	/studies/{study}/series/{series}/instances/{instance}/thumbnail
Frame Thumbnail	/studies/{study}/series/{series}/instances/{instance}/frames/{frames}/thumbnail

The options and restrictions valid for all thumbnail resources are listed in Table 15.

Table 15: Options and restrictions on thumbnail resources

Options	Restrictions
Accept (Request Header Field)	Restricted to image/jpeg, image/gif, image/png, or wildcard representations covering the three rendered media types.
Accept-Charset (Request Header Field)	Not applied as the response payload consists of a bitmap image only.
Accept (Query Parameter)	Restricted to image/jpeg, image/gif, or image/png.
Charset (Query Parameter)	Not applied as the response payload consists of a bitmap image only.
Viewport (Query Parameter)	If viewport width and height are specified, the rendered images will be of the specified dimensions with black background color if padding is necessary. Further parameter values are ignored.
Transfer Syntaxes Supported	The DICOM instance to render as a thumbnail is required to be of a supported Transfer Syntax. For the list of supported Transfer Syntaxes, see Table 16.
SOP Class Restrictions	The DICOM instance to render as a thumbnail needs to be of a supported SOP Class. Rendering is restricted to SOP Classes for the purpose of image data storage - e.g. the DICOM Image Pixel Module exists along with the DICOM Element PixelData (7FE0, 0010).
Image Type Restrictions	Binary images are not supported.
Size Restriction	There are no size limits imposed by the WADO-RS AE. However, depending on system resource and file format restrictions there might be limits to the size of each image.
Frames	The given frame numbers for rendered frames are expected to be in ascending order. If this is not the case, automatic sorting applies.



Note

The first frame of the first instance of the first series is used to render the thumbnail.

C.2.4 WADO-RS supported Transfer Syntaxes

The supported Transfer Syntaxes for all rendered resources and thumbnail resources are listed in Table 16.

Table 16: WADO-RS supported Transfer Syntaxes

Transfer Syntax Name	UID
Implicit VR Little Endian: Default Transfer Syntax for DICOM	1.2.840.10008.1.2
Explicit VR Little Endian	1.2.840.10008.1.2.1
Deflated Explicit VR Little Endian	1.2.840.10008.1.2.1.99
Explicit VR Big Endian	1.2.840.10008.1.2.2
JPEG Baseline (Process 1): Default Transfer Syntax for Lossy JPEG 8-bit Image Compression	1.2.840.10008.1.2.4.50
JPEG Baseline (Processes 2 and 4): Default Transfer Syntax for Lossy JPEG 12-bit Image Compression (Process 4 only)	1.2.840.10008.1.2.4.51
JPEG Lossless, Nonhierarchical (Processes 14)	1.2.840.10008.1.2.4.57
JPEG Lossless, Nonhierarchical, First-Order Prediction (Processes 14 [Selection Value 1]): Default Transfer Syntax for Lossless JPEG Image Compression	1.2.840.10008.1.2.4.70
JPEG-LS Lossless Image Compression	1.2.840.10008.1.2.4.80
JPEG-LS Lossy (Near-Lossless) Image Compression	1.2.840.10008.1.2.4.81
JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90
JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91
RLE Lossless	1.2.840.10008.1.2.5

C.2.5 WADO-RS connection policies

General

All standard RS connection policies apply. There are no extensions for RS options.

C.2.5.1 WADO-RS endpoint URL

`https://{host}/webrenderserver/rs/tenants/{tenant}`

Please replace {host} and {tenant} in the URL with the corresponding host and tenant name.

C.2.5.2 Security

To protect the WADO-RS interface from unauthorized access special attention is paid to security aspects. Only requests over HTTPS are accepted by the WADO-RS AE. On top, requests are required to contain valid authorization information by means of URL signatures (HMAC). Please note that the WADO-RS service is not enabled per default in the *webrenderserver*. For general information about configuration and URL signatures, please consult the *syngo.share* Interface and Integration Manual.

C.2.5.3 Number of connections

The WADO-RS AE itself does not limit the number of simultaneous requests. However, in standard deployments the number of simultaneous HTTP connections might be limited by Traefik.

C.2.5.4 Synchronous and asynchronous requests

The WADO-RS AE only supports synchronous requests.

C.2.5.5 Response status

The response message header contains a HTTP status code indicating success or failure, as listed in Table 17.

Table 17: WADO-RS HTTP status codes

Code	Name	Description
200	OK	The requested resource has been fetched successfully and is returned in the message body.
206	Partial Content	Indicates that for the requested resource (rendered resource) one or more DICOM instances are of a not supported SOP Class and are omitted in the response.
400	Bad Request	The request is malformed. Details are returned in the payload of the response.
403	Forbidden	The WADO-RS AE is refusing action as no or (only) an invalid authentication (URL signature) is present.
404	Not Found	The WADO-RS AE could not find the specified resource. This response status is also used to indicate that no applicable DICOM instance for thumbnail generation was found.
405	Method Not Allowed	Indicates a function parameter keyword other than linear is used in the window query parameter.
406	Not Acceptable	The remote AE did not provide any rendered media type in the accept header field and accept query parameter.
500	Internal Server Error	The WADO-RS AE encountered an error while processing the request.

C.2.5.6 Warning header

The response message header contains a HTTP Warning header field with the number of not rendered instances in case the response status is 206 (Partial Content). Furthermore, if unsupported annotation keywords are given in the request, the unsupported keywords are listed in the HTTP Warning header field.

Distributed by

Siemens Healthineers AG
Siemensstr. 3
91301 Forchheim
Germany

İthalatçı ve Türkiye Yetkili Temsilcisi
Siemens Healthcare Sağlık A.Ş.
Yakacık Caddesi No 111
34870 Kartal İstanbul
Türkiye

Manufacturer

ITH icoserve technology for healthcare GmbH
Innrain 98
6020 Innsbruck
Austria
Phone: +43 512 89059