
Absolute Dose Values in CT

SOMATOM Force

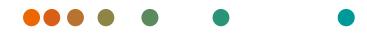

Reference values		Germany ¹	European Union ²	USA ³
Head Routine	CTDIvol[mGy]	60	60	75
Thorax Routine	CTDIvol[mGy]	10	30	N/A
Abdomen Routine	CTDIvol[mGy]	15	35	25

¹Bundesamt für Strahlenschutz (Diagnostische Referenzwerte für diagnostische und interventionelle Röntgenanwendungen vom 22.06.2016).

 $^{^3}$ American College of Radiology (CT Accreditation Program Requirements, Amended 2019).

Default Siemens Protocol	Standard CTDIvol / [mGy]*	Put your current values here
Head Routine	47.32	
Thorax Routine	2.66	
Turbo Flash Lung Low Dose with Tin Filter	0.37	
Abdomen Routine	9.83	

^{*}Values are based on the default protocols of the SOMATOM Force with *syngo* CT VB20 and an average sized patient of 1.75 m and 75 kg. Iterative Reconstruction ADMIRE is used.

In clinical practice, the use of ADMIRE may reduce CT patient dose depending on the clinical task, patient size, anatomical location, and clinical practice. A radiologist and a physicist should be consulted to determine the appropriate dose to obtain diagnostic image quality for the particular clinical task.

²European Guidelines on Quality Criteria for Computed Tomography (http://www.drs.dk/guidelines/ct/quality/htmlindex.htm).

On account of certain regional limitations of sales rights and service availability, we cannot guarantee that all products included in this brochure are available through the sales organization of Siemens Healthineers worldwide.

Availability and packaging may vary by country and is subject to change without prior notice.

Some/All of the features and products described herein may not be available in the United States.

The information in this document contains general technical descriptions of specifications and options as well as standard and optional features which do not always have to be present in individual cases.

The statements by customers of Siemens Healthineers described herein are based on results that were achieved in the customer's unique setting. Since there is no "typical" hospital and many variables exist (e.g., hospital size, case mix, level of IT adoption), there can be no guarantee that other customers will achieve the same results.

The customers cited are employed by an institution that might provide Siemens Healthineers product reference services, R&D collaboration or other relationship for compensation pursuant to a written agreement.

Siemens Healthineers reserves the right to modify the design, packaging, specifications, and options described herein without prior notice. Please contact your local Siemens Healthineers sales representative for the most current information.

Note: Any technical data contained in this document may vary within defined tolerances. Original images always lose a certain amount of detail when reproduced.

The products/features and/or service offerings (here mentioned) are not commercially available in all countries and/or for all modalities. If the services are not marketed in countries due to regulatory or other reasons, the service offering cannot be guaranteed. Please contact your local Siemens Healthineers organization for further details.

Siemens Healthineers Headquaters

Siemens Healthcare GmbH Henkestr. 127 91052 Erlangen, Germany Phone: +49 9131 84-0 siemens-healthineers.com