



Reaching further in X-ray imaging.
Mobile. Digital. Wireless.

Mobilett Mira

Answers for life.

SIEMENS

How can I always stay at the cutting edge of mobile X-ray technology?

How can I get stationary-like image quality, the mobile way?

How can I combine mobility with easy handling, first-class images, and the utmost connectivity?

Mobilett Mira

Reaching further than ever before in X-ray imaging

The name Mobilett stands for mobility: the system is easy to maneuver, fast to use, powerful, and flexible in every situation. Today, however, just being mobile isn't enough.

That's why we take digital mobile X-ray imaging to the next level – with our new generation Mobilett Mira. The fully integrated system provides the answers to your everyday challenges.

For you, that means mobility combined with:

- simple operation
- intelligent networking
- uncompromised image quality
- a high level of security and availability.

Experience wireless freedom and flexibility in the mobile X-ray sector. Experience Mobilett Mira.

What makes the new really new?

With the new Mobilett Mira, digital mobility goes wireless – combining increased freedom of movement, simple handling, and excellent image quality.

Get to know the innovations of Mobilett Mira and experience digital imaging at its best!



Wireless detector

Wireless and more – from image quality to mobile handling.

Whether in daily handling or in image transfer, Mobilett Mira gives you wireless flexibility. The result: safer operation, faster digital data transfer, and no cables to trip over, or become entangled. The image quality is comparable to that of high-end stationary systems, with a resolution of seven million pixel for confident diagnoses.

You can take advantage of all these benefits around the clock. The batteries of both the detector and the entire system are automatically charged with a single central charging mechanism. In fact, everything is as simple and straightforward as you need it during your mobile routine.

At a glance

- Seamless integration of wireless technology into digital imaging
- Outstanding detector resolution
- Complete imaging freedom in your daily routine



Reaching further – wireless, everywhere you go

Rotating swivel arm

Flexibility and cleanability – in a single arm system.

The rotating swivel arm with its industry-leading arm range gives you the flexibility you’ve always wanted for a mobile X-ray system. The counterbalanced design lets you examine every region of interest in every position and from almost every angle, even with patients who can’t move. With the combination of industry-leading arm range and the new 90 degree freedom of lateral arm rotation to the left and right, a wide range of examination positions is possible.


In short: simple and fast positioning for any patient environment. The fully integrated cable design ensures greater sterility, and the smooth surfaces are very easy to clean and disinfect. This is one arm that puts all the advantages in your hand.

At a glance

- Industry-leading arm range
- New degree of freedom by means of lateral arm rotation
- Hygienic design due to integrated cables



Reaching further – more flexibility, everywhere you go



Highest imaging power per footprint

Uncompromised image quality – wherever you go.

With Mobilett Mira you get the highest imaging power in the smallest space. The high output power of up to 35 kW or 450 mA gives you remarkably sharp images without motion blurring, images that are fully comparable to those of high-end stationary systems – yet with extremely short exposure times of as low as one millisecond, all in a very compact unit design.

The state-of-the-art WLAN detector with high-sensitivity CsI as a scintillator material provides first-class image quality under any conditions. That means lower radiation dose for your patients and even more diagnostic certainty for you.

UPTIME Services

Peace of mind – reliable running.

During the hectic day-to-day hospital routine, you need a system that is reliable and always up-to-date. A workhorse like the Mobilett Mira: tough, durable, ready for anything. And not just in terms of daily use, but also in terms of data security. Regular anti-virus updates give you solid protection against malware attacks and ensure you the uptime you need.

Security and availability are indispensable for a mobile X-ray unit. With Siemens Remote Service (SRS), you have a reliable, proactive partner. You take care of your patients, while we ensure uninterrupted system availability, online with remote diagnosis, and on-site at your facility. Thanks to remote update handling, you also benefit from the latest Siemens technology developments. Mobilett Mira brings together dependability with advanced system design.

At a glance

- High imaging output for stationary-like images
- Excellent image quality in a compact system
- High-resolution FD for low dose applications



At a glance

- Virus protection
- Remote Diagnosis
- Remote Update Handling




Reaching further – brilliant images, everywhere you go

Reaching further – fully supported, everywhere you go.

How does Mobilett Mira improve **your workflow?**

Mobilett Mira combines digital mobility with wireless technology, making your daily work with patients easy and straightforward. Since it's always ready to go, you can reach more patients and overcome nearly all spatial barriers. Thanks to its friendly giraffe design, it also creates a relaxed patient environment.

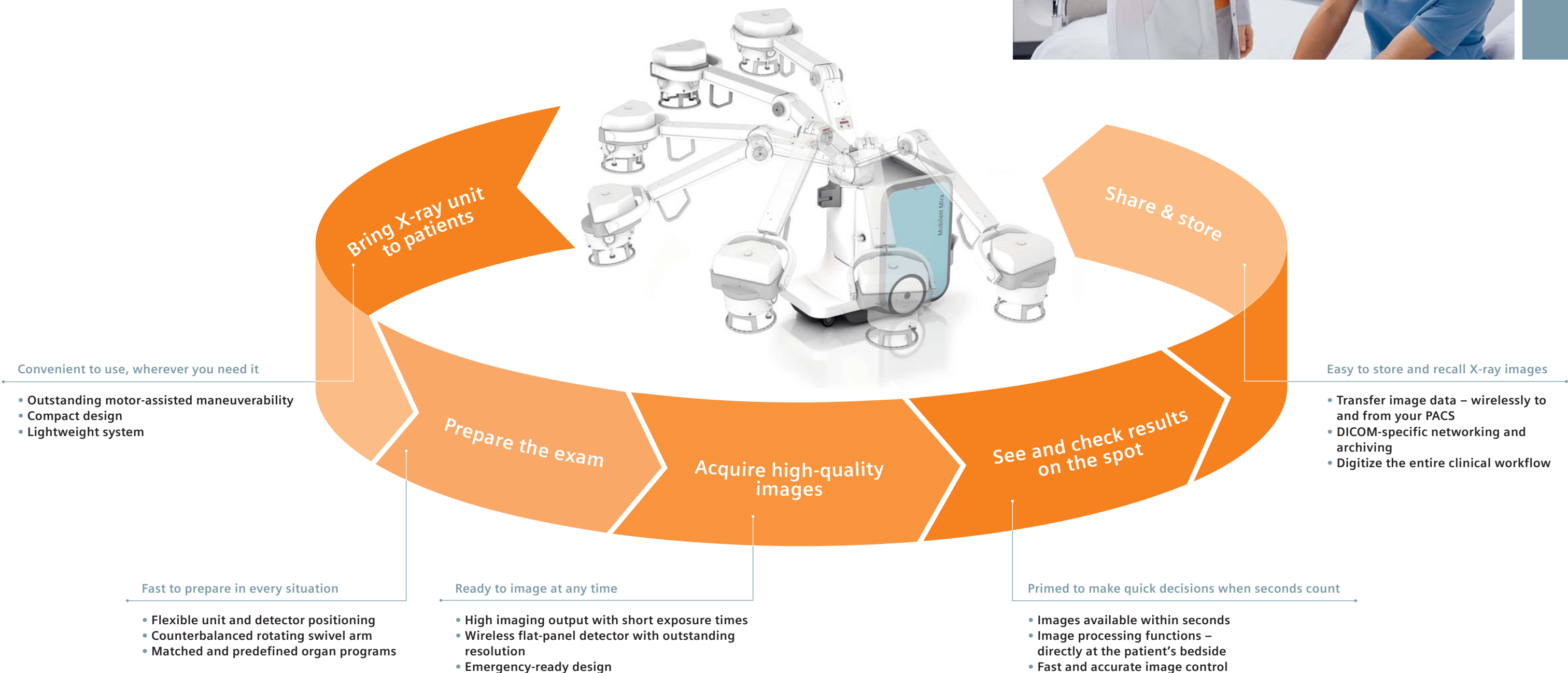

Mobile X-ray imaging at its best



Mobilett Mira helps to make your day simpler and your diagnoses even more certain. The system is the result of many years of experience in mobile X-ray imaging, during which time we have consistently developed the Mobilett family of products.

Wireless freedom, outstanding flexibility, and fully integrated digital mobility make this system unique. Imaging is also more comfortable for patients, the friendly giraffe design creates a relaxed atmosphere and makes the exam seem less threatening – for children and adults alike.

Not only have we advanced the technology of X-ray imaging, we've also made it more environmentally friendly. In addition to being almost free of hazardous internal components, Mobilett Mira has a battery-saving ECO-Mode. This lets you reduce your environmental footprint and use up to 25 percent less energy.

Mobilett Mira is your new mobile X-ray system – reaching further than ever before in X-ray imaging.



This product does not yet have a certificate of conformity. This product is not commercially available in the United States (U.S.A.).

The information in this document contains general descriptions of the technical options available and may not always apply in individual cases.

The required features should therefore be specified in each individual case at completion of contract. Siemens reserves the right to modify the design and specifications contained herein without prior notice. Please contact your local Siemens sales representative for the most current information.

Original images always lose a certain amount of detail when reproduced.

In the interest of complying with legal requirements concerning the environmental compatibility of our products (protection of natural resources, waste conservation), we recycle certain components. Using the same extensive quality assurance measures as for new components, we guarantee the quality of these recycled components.

Please find fitting accessories:
www.siemens.com/medical-accessories

Local Contact Information

In Japan

Siemens-Asahi
Medical Technologies Ltd.
Takanawa Park Tower 14F
3-20-14, Higashi-gotanda
Shinagawa-ku
Tokyo 141-8644
Phone: +81 3 5423 8510

In Asia

Siemens Pte Ltd
The Siemens Center
60 MacPherson Road
Singapore 348615
Phone: +65 6490 8182

In Germany

Siemens AG
Healthcare Sector
Karlheinz-Kaske-Str. 2
91052 Erlangen
Phone: +49 9131 84-0

In the USA

Siemens Medical Solutions USA, Inc.
51 Valley Stream Parkway
Malvern, PA 19355
USA
Phone: +1 610 448 4500

Global Business Unit

Siemens AG
Medical Solutions
X-Ray Products
Henkestr. 127
91052 Erlangen
Germany
Phone: +49 9131 84-0

Global Siemens Headquarters

Siemens AG
Wittelsbacherplatz 2
80333 Muenchen
Germany

Global Siemens Healthcare Headquarters

Siemens AG
Healthcare Sector
Henkestr. 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
www.siemens.com/healthcare

Legal Manufacturer

Siemens AG
Wittelsbacherplatz 2
DE-80333 Muenchen
Germany

www.siemens.com/mira