



Aina

Blood Monitoring System

Enhance your clinic practice with the Aina Station portable lab system. Perform a full panel of metabolic tests while your patients wait.



HbA1c



Glucose



Lipids Profile



Hemoglobin

Marketed by
Siemens Healthcare Private Limited
Unit No. 9A, 9th Floor, North Tower, Godrej One, Pirojshanagar,
Easter Express Highway, Vikhroli East, Mumbai - 400 079
Tel.: +91 22 33700600 | Email: hc_contact.india@siemens-healthineers.com
Helpline for business queries: 1800-209-1800
CIN: U74999MH2015PTC264859


SIEMENS
Healthineers



Aina

Blood Monitoring System

- Real-time, lab-quality results within minutes
- Integrated patient registration and check-in
- Connected cloud-based EMR and doctor app

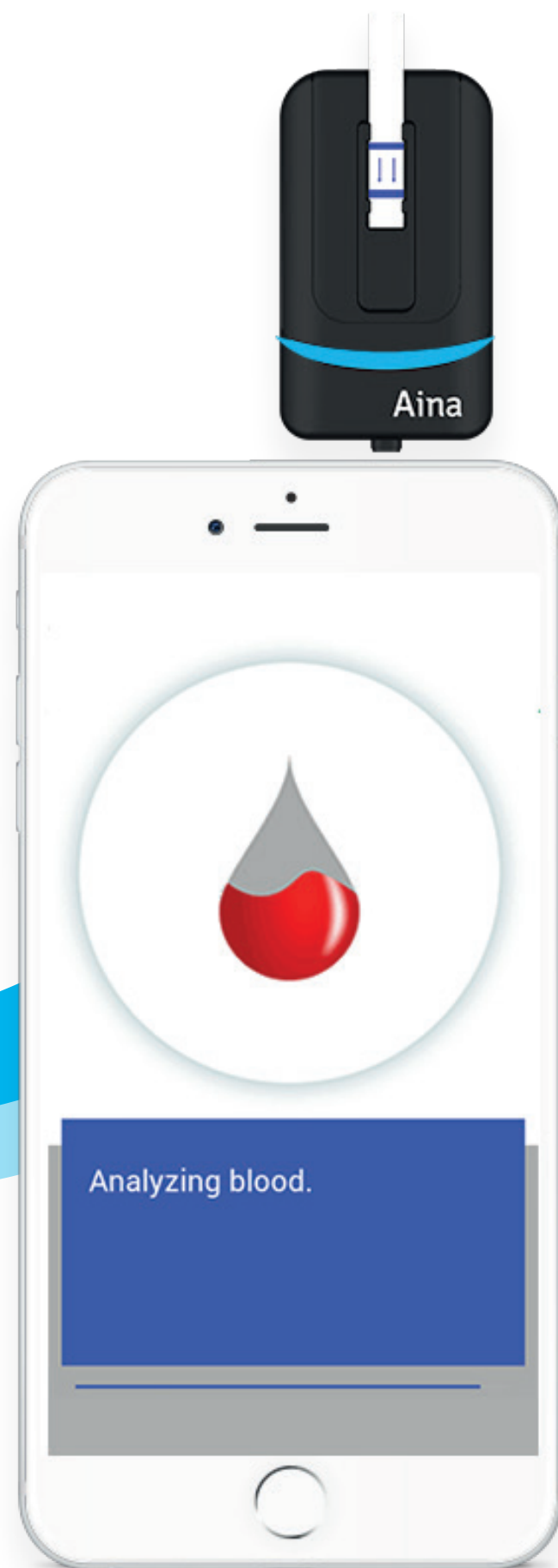

Broad and expanding range
of tests with **one platform**

	HbA1c	Glucose	Total Cholesterol	HDL Cholesterol	Triglycerides	Hemoglobin
Measuring Range	4-15%	10-500 mg/dL	100-300 mg/dL	25-85 mg/dL	50-450 mg/dL	7-23 g/dL
Test Time	3 mins	Approx. 5 secs	Approx. 2 min	Approx. 2 min	Approx. 2 min	30 secs
Sample Size	5 µL	2 µL	15 µL	15 µL	15 µL	10 µL
Ref. Method	Variant II Turbo (Bio-Rad)	2300D STAT Plus (YSI)	Dimension RxL Max (Siemens)	Dimension RxL Max (Siemens)	Dimension RxL Max (Siemens)	LH750 (Beckman Coulter)
Clinical Accuracy	97.6% samples within 10% bias	95.7% samples within 15% bias	100% samples within 20% bias	99% samples within 15 mg/dL bias	95% samples within 15% bias	97.6% samples within 20% bias
R ²	0.972	0.9941	0.973	0.969	0.970	0.964

Marketed by

Jana Med Tech Pvt. Ltd
Bangalore Bioinnovation Center, Bangalore Helix
Biotech Park, Electronic City, Phase I,
Bangalore - 560100, India
www.janacare.com