


ACUSON Juniper ultrasound system

Everything you need, nothing you don't

siemens-healthineers.com/juniper

*The clinical overlay is not that of the individual pictured.
It was modified for better visualization.*

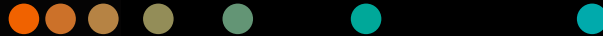

SIEMENS
Healthineers

Helping you transform care delivery

Higher patient volumes. Larger patient sizes. Limited diagnostic certainty. In the face of these trending challenges, you need a reliable ultrasound system that can consistently deliver high-quality images across disease states – a system with optimized accessibility, expanded applications, and ease of use.


Created with these needs in mind, the ACUSON Juniper empowers you to image every patient, knowing you have the clinical information necessary for confident decision-making. Small on the outside, yet big on the inside, the system is designed to adapt to your everyday clinical and workflow challenges.

When form and function meet, the result is the ACUSON Juniper, an ultrasound system that is smaller, more powerful, and capable of accommodating virtually all patients.




The ever growing burden of liver disease

Worldwide 1.9 billion adults are overweight¹ and 600 million are obese¹

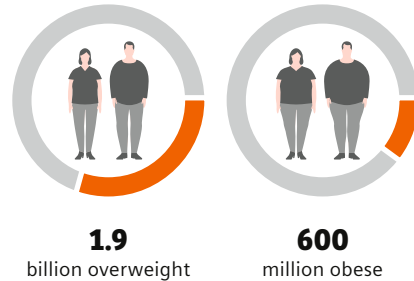

Ultrasound professionals are faced with a patient population that is increasingly harder to scan. Imaging at greater depths is needed to provide higher levels of diagnostic confidence.


Globally, liver cirrhosis has increased by 20 percent in the last 10 years¹

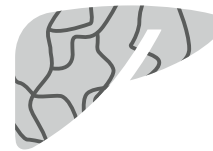

20%
increase of cirrhosis

Millions of people are at risk of developing liver disease from obesity, alcohol abuse, and Hepatitis. However, patients usually show no symptoms until the very late stages. Using ultrasound as a non-invasive test for liver disease can help assess liver tissue stiffness within minutes.

Liver biopsies are costly and carry potential complications

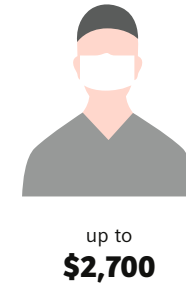

Non-invasive technologies, such as shear wave elastography, produce instant results and potentially reduce healthcare costs by avoiding possible complications from infections caused by interventional procedures².

Everything you need, nothing you don't

Small on the outside, yet big on the inside


ACUSON Juniper is powered by a completely new platform, designed from the ground up to offer high-fidelity transmit and receive acoustic signals that greatly reduces noise and offers premium image quality and industry leading elasticity solutions.

Where form and function meet

Performance comes at no cost to the functional design of the system. ACUSON Juniper offers one of the smallest market footprints in its class, effortless maneuverability, 5 active transducer ports and a simple, intuitive and highly customizable control panel.



Where form meets function

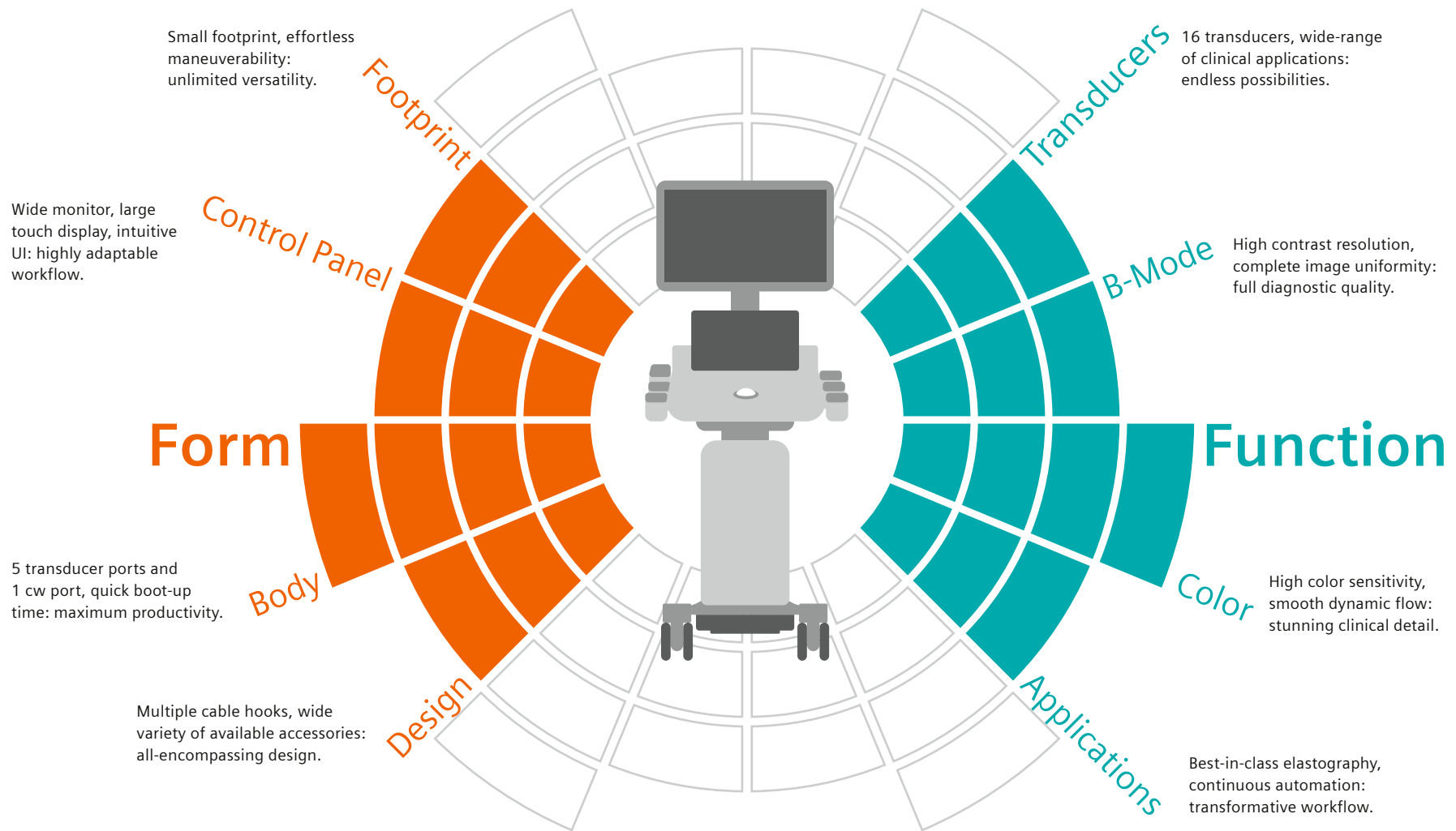

Focus on what matters most

Move easily between exams, patients, and spaces while still benefiting from remarkable system performance with the ACUSON Juniper system. The 6 ports, including 5 active transducer ports and 1 CW port, support multiple clinical applications as you move about your day. Designed with intent, the

system adjusts easily to individual working preferences – pull the lightweight system from the front or back as you move between exams, tilt the monitor up or down as needed for your height, or rotate the control panel left or right depending on the patient bed location.



Control Panel

Support greater workflow personalization with a simple, intuitive control panel and tiltable, highly customizable 13.3-inch touch screen. Able to be optimized for space availability, convenience, and user preference, the system's 90° right and left rotation allows flexibility and improved transducer connector positioning.



The Footprint

ACUSON Juniper offers powerful capabilities in one of the industry's smallest ultrasound systems. Able to fit in virtually any hospital room, the system is up to 36% smaller and weighs an average of 27% less than other systems in its class.


ACUSON Juniper



Vendor A
+14% larger



Vendor B
+31% larger

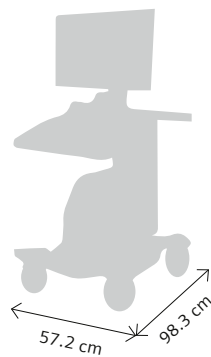

Vendor C
+36% larger



Up to **36% smaller** than other systems in the same class, the **ACUSON Juniper** fits into virtually any hospital room. Weighing **only 76 kg (168 pounds)**, making it one of the lightest and easiest to maneuver Ultrasound system in its class.



21.5-inch LED monitor highly adjustable transition for optimal ergonomic positioning toward, away and side-to-side.



13.3-inch high resolution tiltable touch panel with rapid response technology, allows for quick and accurate key selection.



90° right and left rotation allows the user the flexibility of choosing the transducer connectors position.




Front and back handle for effortless maneuverability.



Individual Castor lock control on all four wheels.

Focus on what matters most

Among the quietest ultrasound systems in its class on the market, the ACUSON Juniper allows a comfortable room environment and redefines adaptability, meeting staff preferences with ease so you can focus on what matters most – your patient.

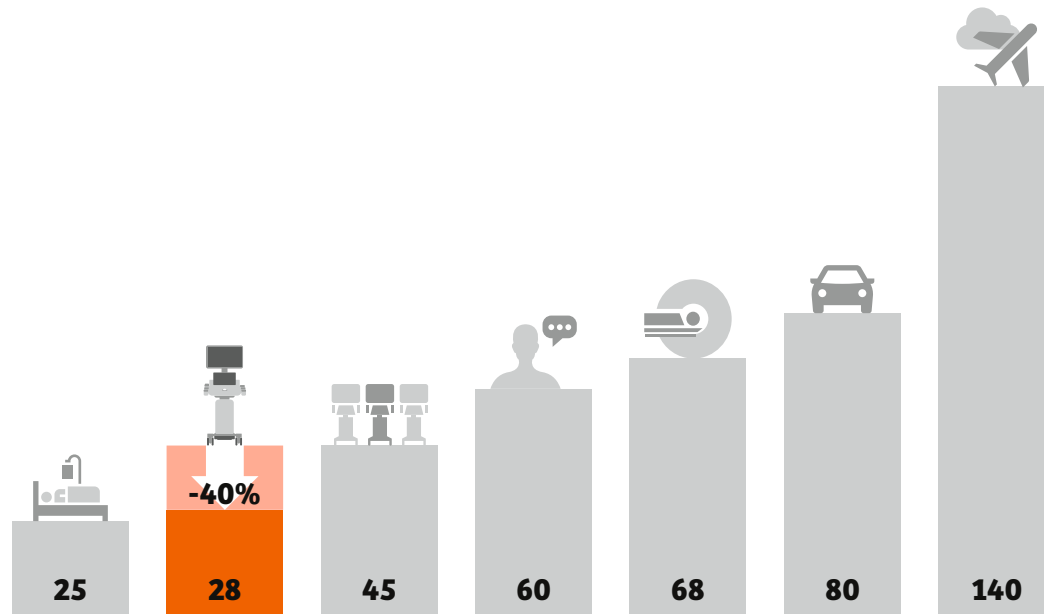

The sound of silence. ACUSON Juniper is **40% quieter** than the average ultrasound system. Sound pressure in decibels (dB)



Design

The control panel supports a clean and sterile environment with its customized disposable skin. Several cable hooks are available to support better cable management and efficiency.



Body

With 5 active transducer ports and 1 pencil port, the ACUSON Juniper is always ready to scan with little to no setup time. This system is designed for a stress-free environment with a sound pressure of only 28 decibels.



Customized disposable skin keeps the control panel clean and sterile.



Purposeful cable management hooks – several hooks available for a clean, cable free environment.



Everything made with intent, including utilities storage basket, magnetic storage bin, and transducer holders.



Efficient workflow with 5 active transducer ports and 1 CW port.



ECG cable connector, DVD writer and physio port.

Image every patient


Routinely delivering high-quality care means scanning virtually every patient regardless of their size, weight, or condition. With its clinical versatility and industry-leading tissue elasticity imaging, the ACUSON Juniper enables you to expand your ultrasound clinical and service offerings across virtually all patient and case types.

The system's wide variety of clinical solutions and fast image acquisition can improve diagnostic confidence across all departments, from radiology, interventional radiology, and urology to cardiovascular, orthopedic, and obstetrics/gynecological imaging.



Transducers

Cover a wide range of clinical needs with 16 transducers that deliver the versatility to expand clinical offerings.

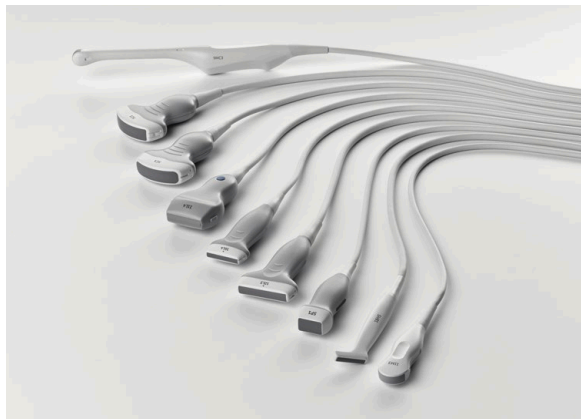

Shared Services Applications



Women's Health Applications




Cardiac Applications




B-Mode

High contrast resolution, complete image uniformity: full diagnostic quality.

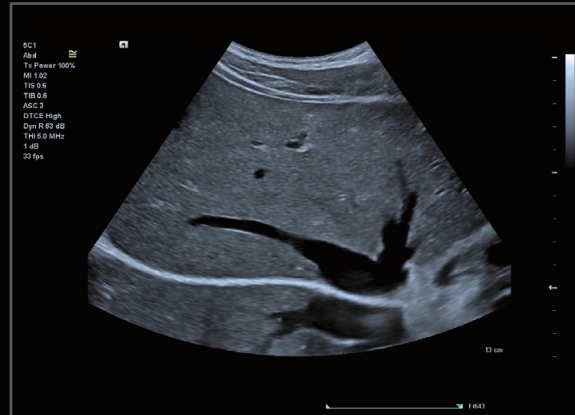


Complete homogeneous liver texture from the near to the far field with a well-balanced pixel pattern and sharp border delineation.

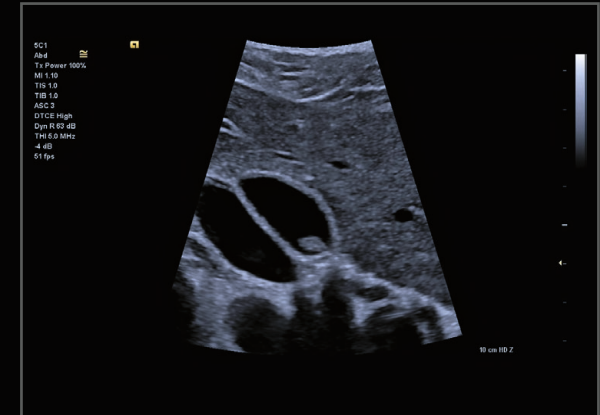


Improved tissue differentiation and border detection, as with this gall bladder mass, using Dynamic TCE tissue contrast enhancement technology.



This image of a fetal profile shows excellent contrast resolution and visualization of subtle tissue differences using eSiImage Multiparametric Optimization, without having to adjust any of the gain controls.

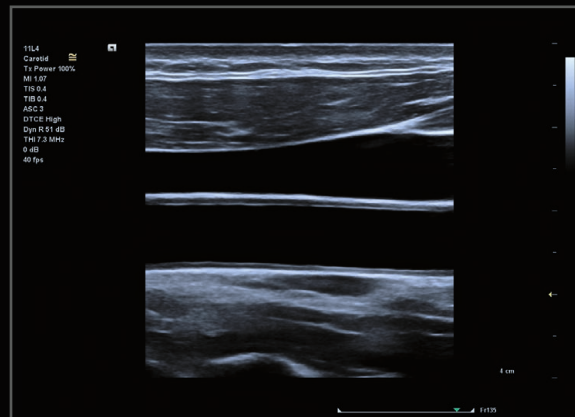

A simple and intuitive control panel allows the users to easily adjust the imaging parameters, resulting in increased speed and work flow efficiency.

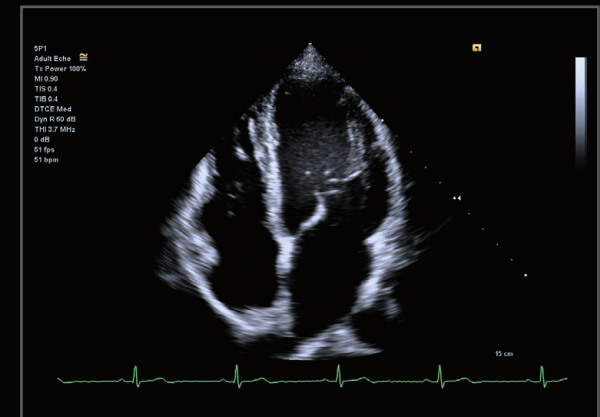


Single crystal technology improves diagnostic confidence while performing dynamic cardiac scans as demonstrated in this view of the four chamber heart.

Image every patient

The ACUSON Juniper delivers consistency and repeatability for your clinical environment by guiding the user on what views are required and what else they need to do to get a complete study. Reduce motion

artifacts and noise while simultaneously enhancing color sensitivity with Dynamic Persistence and color Auto Flash Artifact Suppression technologies.

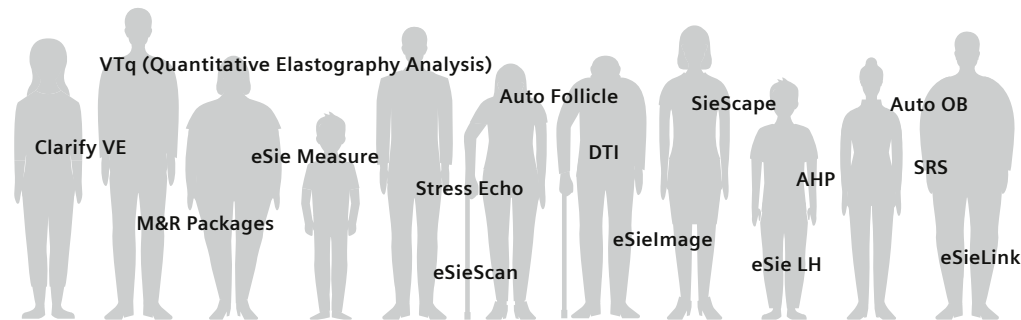

ACUSON Juniper Advanced applications for maximum productivity and confidence.




Color

Improve diagnostic confidence with exceptional color resolution and with improved physiological representation of color while reducing color flash artifact.




Applications


Offers the ability to measure tissue elasticity imaging with a click of a button using our best-in-class Virtual Touch quantification and eSie Touch elastography. The system helps to reduce operator variability, improves accuracy and gives greater consistency in reporting through several automated tools.

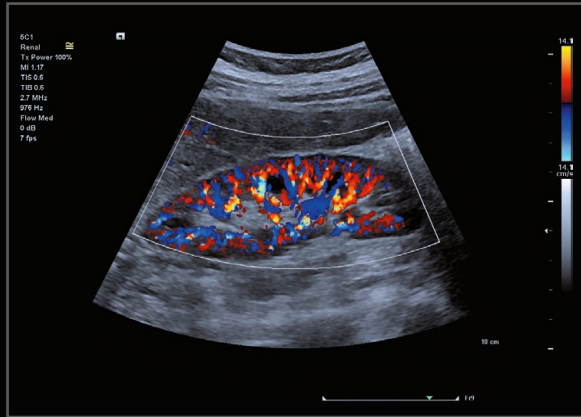


The perfusion of the renal cortex can be seen here using our color Doppler with Auto Flash Artifact Suppression technology.

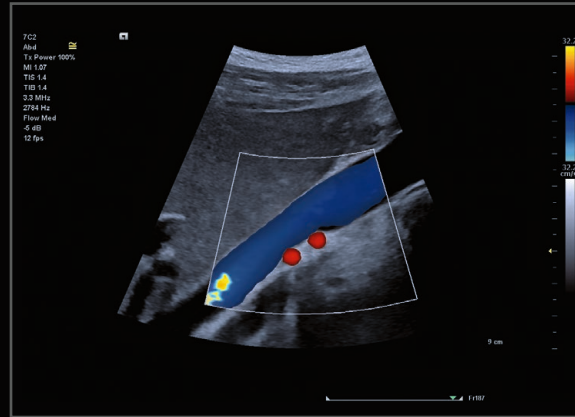


Dynamic Persistence and Auto Flash Artifact Suppression technologies help reduce tissue motion and color Doppler artifact, shown in this image of the IVC and duplicated right renal artery.

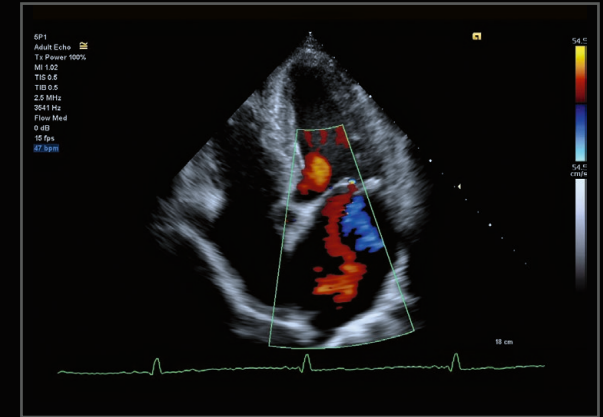


A small mitral valve regurgitation is seen with our True Fidelity color Doppler technology.

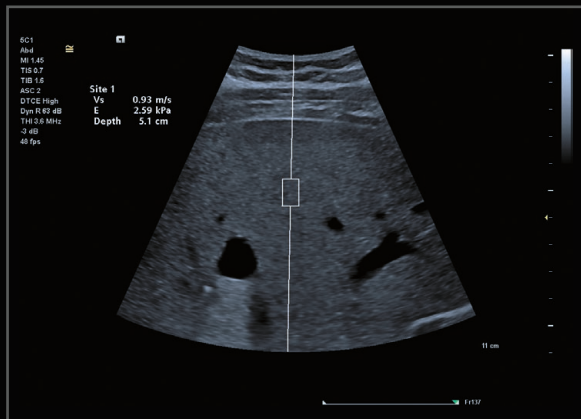

Reproducible, reliable and detailed tissue stiffness information can be quickly and easily obtained using our one touch Virtual Touch quantification technology.

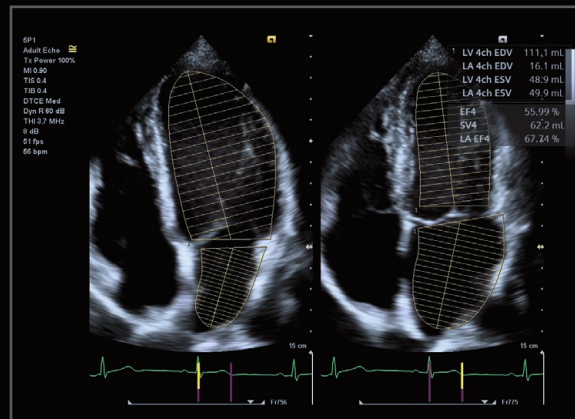

Automated measurement of volume and ejection fraction for LV and LA function in 2D echo using eSie Left Heart.

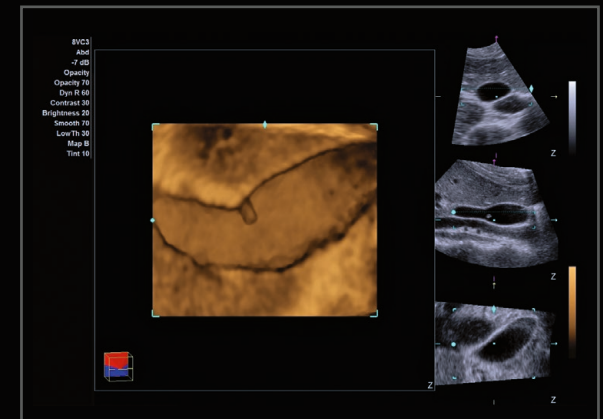

3-Scape 3D shows a detailed representation of the gallbladder pathology seen here.

Where form meets function

Functional Performance

- Shear wave elastography imaging
- Comprehensive measurements package
- 16 new transducers to cover all clinical applications
- New front-end, beamformer and back-end engines for improved contrast resolution, plunkability, sensitivity and specificity

13.3" touch display

Smart UI for quicker scan settings and protocols

Intuitive control panel

Customizable height, angle and keys

Quiet, stress-free environment

Bedroom-like sound level, 40% quieter than average

Every detail designed to maximize productivity

Multiple storage bins available as well as space for easy installation of peripherals



21.5" LED monitor

Highly adjustable for maximum operator adaptability

Lightweight maneuverability

One of the easiest systems to maneuver, weighing only 76 kilos (168 pounds), and with both rear and front handles

Purposeful cable management hooks


Several hooks available for a clean, cable free environment

Multiple transducer ports

Efficient workflow with 5 active transducer ports, 1 pencil cardiac port and 2 physio ports

Reliable support

eSieLink:
High uptime
and increased
user support

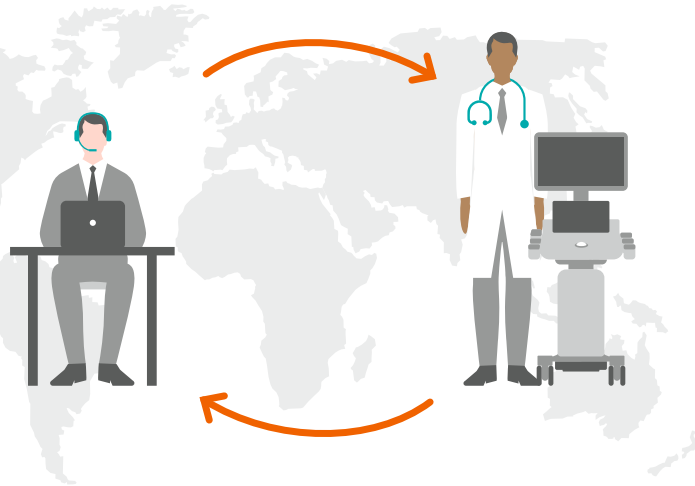

Service: Remote assistance technology

You need to know that you have access to support when and where you need it. With eSieLink Remote Assistance technology, you get just that so you can work more confidently. eSieLink enables fast problem resolution and overall increased productivity.

Cloud-based network for imaging fleet management

teamply is a cloud-based network that brings together healthcare professionals in order to advance medicine and human health as a team effort. Be it an institution, hospital chain, or integrated delivery network, transparency of performance is key to its success. teamply brings this transparency to the imaging fleet. With its built-in multi-vendor support, it grants instant access³ to fleet statistics. And more crucially, it empowers healthcare professionals to identify improvement potential on all levels of execution.



The products/features mentioned in this document may not be commercially available in all countries. Due to regulatory reasons their future availability cannot be guaranteed. Please contact your local Siemens organization for further details.

Standalone clinical images may have been cropped to better visualize pathology.

ACUSON Juniper, eSieLink, eSie Touch, eSie Left Heart (LH), eSieScan, eSie Measure, Clarify VE, eSieImage, 3-Scape, and Virtual Touch quantification (VTq) are trademarks of Siemens Medical Solutions USA, Inc.

At Siemens Healthineers, our purpose is to enable healthcare providers to increase value by empowering them on their journey towards expanding precision medicine, transforming care delivery, and improving patient experience, all enabled by digitalizing healthcare.

An estimated 5 million patients globally everyday benefit from our innovative technologies and services in the areas of diagnostic and therapeutic imaging, laboratory diagnostics and molecular medicine, as well as digital health and enterprise services.

We are a leading medical technology company with over 170 years of experience and 18,000 patents globally. With more than 48,000 dedicated colleagues in 75 countries, we will continue to innovate and shape the future of healthcare.

Endnotes

¹ World Health Organization (WHO), 2010

² Franciscus, Alan, 2014. Hepatitis C Support Project, hcvadvocate.org/hepatitis/factsheets_pdf/biopsy.pdf

³ Prerequisites include: wireless connection to clinical network, meeting recommended minimum hardware requirements, and adherence to local privacy and security regulations.

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen, Germany
Phone: +49 9131 84-0
siemens-healthineers.com

Legal Manufacturer

Siemens Medical Solutions USA, Inc.
Ultrasound
685 E. Middlefield Road
Mountain View, CA 94043, USA
Phone: 1-888-826-9702
siemens-healthineers.com/ultrasound