


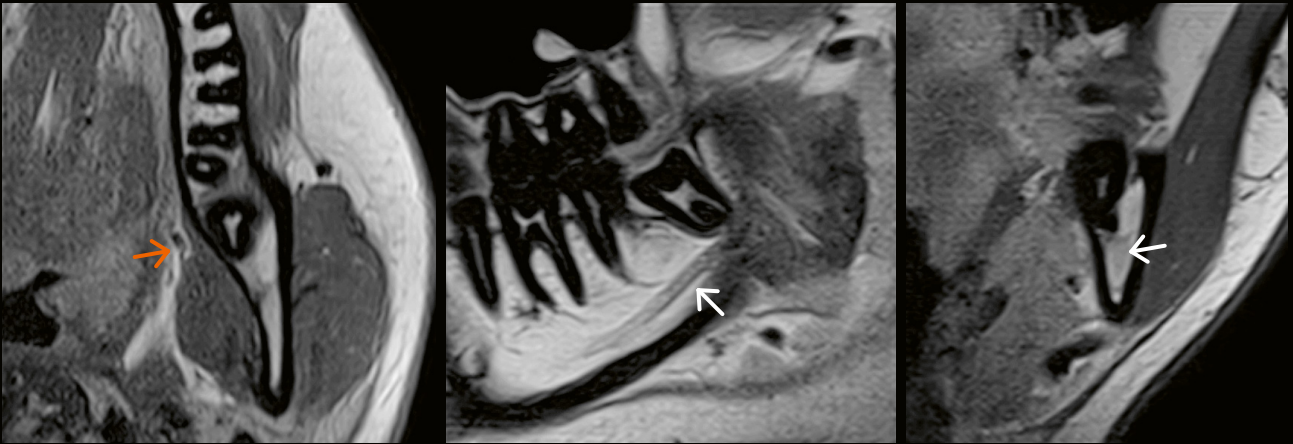

The Journey of Dental-Dedicated MRI at a School of Dentistry: The Aarhus Experience

Jennifer Christensen and Rubens Spin-Neto

(Department of Dentistry and Oral Health, Aarhus University, Aarhus, Denmark)

MRI for the Dental Practice. A Technology Finds its Way

Miriam Keil

(Siemens Healthineers, Erlangen, Germany)

The Journey of Dental-Dedicated MRI at a School of Dentistry: The Aarhus Experience

Jennifer Christensen and Rubens Spin-Neto

Department of Dentistry and Oral Health, Aarhus University, Aarhus, Denmark

Abstract

Dental-dedicated magnetic resonance imaging (ddMRI) has recently emerged as a promising radiation-free imaging modality for dentomaxillofacial diagnostics. By combining a lower-field MRI system, a dedicated radiofrequency coil, and optimized pulse sequences, ddMRI enables simultaneous visualization of dental structures, bone, soft tissues, and inflammatory processes. Aarhus University was among the first institutions worldwide to implement and systematically investigate this technology within its school of dentistry. Through a series of experimental *ex vivo* and clinical studies, the researchers explored sequence optimization, workflow integration, and diagnostic feasibility across multiple dental disciplines. This perspective article summarizes the first years of ddMRI implementation in Aarhus, highlighting technological development, patient workflow, sequence optimization, and emerging clinical applications. The Aarhus experience illustrates both the potential and the current limitations of ddMRI and outlines key research that will be necessary for integrating MRI into routine dentomaxillofacial imaging.

Introduction

Diagnostic imaging is fundamental to modern dentistry and plays a key role in diagnosis, treatment planning, and follow-up across dental specialties [1, 2]. Conventional dental imaging relies predominantly on radiographic techniques, including intraoral radiography, panoramic imaging, and cone beam computed tomography (CBCT) [2]. These modalities provide excellent spatial resolution for mineralized tissues. However, they also expose patients to ionizing radiation and provide only limited information about soft tissues and inflammatory processes [1–3].

Magnetic resonance imaging (MRI) offers radiation-free imaging and excellent soft-tissue contrast. It has

therefore long been considered a potentially valuable modality for dental imaging [4, 5]. However, MRI has historically played a limited role in dentistry due to infrastructural requirements, cost, and the technical challenges of imaging highly mineralized dental tissues [4, 5].

Recent technological developments have resulted in the emergence of dental-dedicated MRI (ddMRI), which combines lower-field MRI scanners (i.e., 0.55T), dedicated dental surface coils, and optimized pulse sequences tailored for dentomaxillofacial imaging [1–3]. These innovations allow improved visualization of dental structures and surrounding tissues while maintaining the intrinsic advantages of MRI, including multiplanar imaging and radiation-free acquisition [1–3].

The Section for Oral Radiology and Endodontics in the Department of Dentistry and Oral Health at Aarhus University, Denmark, was among the first clinical and academic environments to systematically explore this technology. Over the past several years, a series of technical *in vitro*, *ex vivo*, and clinical studies have evaluated ddMRI with respect to sequence development, workflow integration, and diagnostic feasibility. This perspective article summarizes the early “Aarhus experience” and discusses how ddMRI might influence the future of dentomaxillofacial imaging.

Technological foundations of dental-dedicated MRI

Dental-dedicated MRI differs from conventional MRI systems in several important aspects [3]. First, the only ddMRI system available operates at a lower magnetic field strength (0.55T). This system has a smaller footprint, requires less structural preparation of the building, and consumes fewer resources such as electricity and water. Also, lower field strengths have the major advantage of reducing the magnitude of susceptibility artifacts caused

by metallic restorative materials (e.g., fillings and implants) and air-tissue interfaces, which are frequently encountered in the oral cavity [6].

Second, ddMRI employs dedicated surface coils specifically designed to fit and cover the dentomaxillofacial region. These coils maximize signal-to-noise ratio in the region of interest and allow acquisition of high-resolution datasets approaching the spatial resolution required for dental diagnostics, which is typically higher than the image resolution used for “regular” MRI (e.g., 0.2–0.4 mm, versus the usual 0.5–1.0 mm) [7].

Third, pulse sequences must be optimized for dental tissues. Teeth and cortical bone exhibit low proton density and short relaxation times, requiring tailored acquisition parameters. Proton-density-weighted sequences have proven useful for anatomical depiction, whereas fat-suppressed sequences such as short-tau inversion recovery (STIR) highlight inflammatory processes such as edema associated with periodontal and periapical pathology.

Together, these developments enable ddMRI to visualize dentomaxillofacial anatomy with increasing spatial resolution while maintaining the key advantage of radiation-free imaging.


Implementation of dental-dedicated MRI in Aarhus

The introduction of ddMRI at Aarhus University was the earliest implementation of the technology within a dental academic setting. The system used in the Aarhus studies

consisted of a lower-field, 0.55T MRI scanner equipped with a dedicated seven-channel surface dental receive coil (MAGNETOM Free.Max Dental Edition ddMRI system; Siemens Healthineers, Erlangen, Germany, and RAPID Biomedical, Rimpf, Germany).

The implementation required adapting conventional dental imaging workflows — including patient screening procedures and safety protocols — to MRI-specific requirements. It also involved developing optimized acquisition protocols that reflected the clinical demand (i.e., the diagnostic questions) and were tailored to the dental practice (e.g., point-of-care availability, ease of operation, shorter acquisition times). The process started in late 2022 and ended successfully with the finalization of the clinical trials that provided the evidence necessary for the system to achieve CE and FDA clearance in 2025.

Typical examinations begin with a scout acquisition followed by targeted sequences optimized for anatomical visualization, inflammatory assessment, and tissue differentiation. Acquisition times typically range from 10 to 25 minutes, depending on the clinical indication. These protocols were refined through a combination of phantom experiments, *ex vivo* studies, and clinical pilot investigations (Fig. 1). This translational workflow allowed rapid feedback between technical development and clinical feasibility. The current workflows reflect the basic steps for implementing the new diagnostic modality; they allow room for improvement, mostly by making the protocols more targeted to each application and avoiding redundancy.

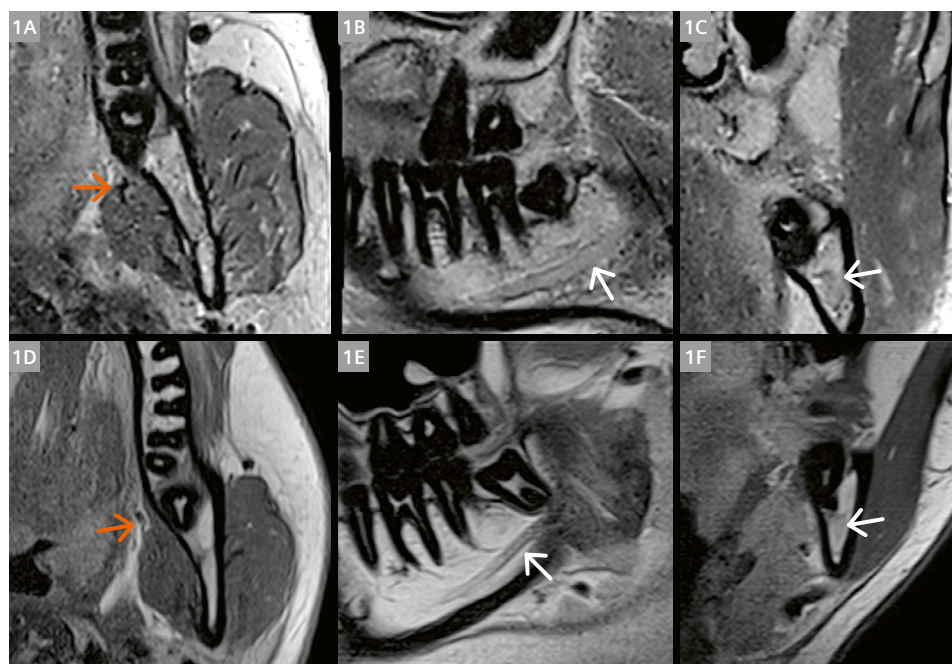

1 Proton-density-weighted images of impacted inferior third molars. The orange arrows show the submandibular ganglion and its connection to the lingual nerve; the white arrows show the mandibular canal (1A and 1D, axial view; 1B and 1E, sagittal view; 1C and 1F, coronal view). The images in the top row reflect the early pulse sequence development and had an acquisition time of 5 minutes and 12 seconds each. The images in the bottom row reflect the current status and had an acquisition time of 2 minutes and 38 seconds each.



2 Proton-density-weighted images of the TMJ region. The orange arrows show the articular disc. The white arrows show the articular space. **(2A)** shows the disc anteriorly displaced while the patient had a closed mouth. **(2B)** shows the disc "recaptured" at the regular position when the patient opened their mouth.

Patient experience and clinical workflow


Introducing ddMRI into dental diagnostics requires consideration of patient acceptance and workflow integration. In ddMRI examinations, patients are positioned supine in the scanner, and the dental surface coil is positioned around the jaw. As with any form of dental-related imaging, the patients were not asked to change clothes.

In the workflow established at Aarhus University, patients were screened for MRI contraindications, including pacemakers, ferromagnetic implants, and severe claustrophobia. Early clinical pilot studies demonstrated that ddMRI examinations could be integrated into dental imaging workflows and that patient tolerance was generally good [3, 8–12].

Although scan times are longer than those of radiographic imaging, early experience suggests that ddMRI may be feasible for selected clinical indications where radiation-free imaging provides additional diagnostic value. The proposed workflows cover indications that commonly lead to referrals for diagnostic imaging in dental patients: temporomandibular joint (TMJ) assessment (Fig. 2), assessment of endodontic and periodontal diseases (Fig. 3), planning for the extraction of third molars (and other teeth), and orthodontic treatment planning [3]. Furthermore, basic oral anatomy and inflammation visualization are included in the proposed workflows.

Sequence development and imaging optimization

Fine-tuning and optimizing sequences has been a central focus of the research so far. Proton-density turbo spin-echo (PD-TSE) sequences have been widely used for anatomical depiction of dental structures. Three-dimensional proton-density sequences (PD-SPACE) enable isotropic image

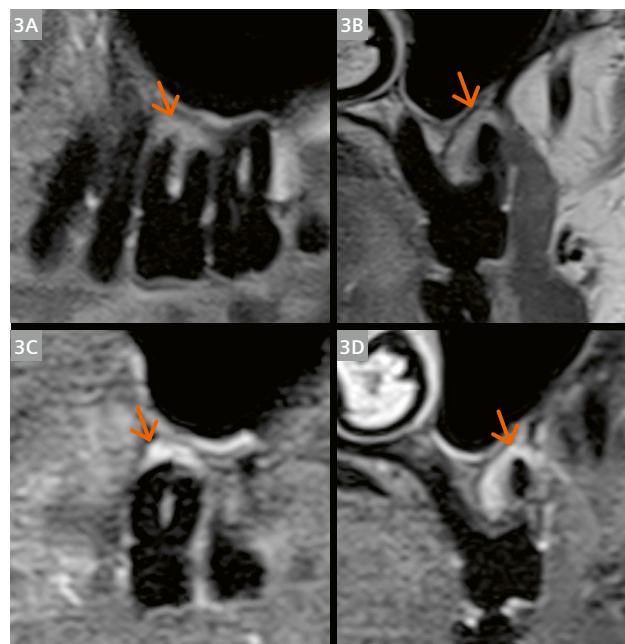

3 Proton-density-weighted images without fat suppression (**3A** and **3B**) and with fat suppression (**3C** and **3D**) of an upper molar with apical periodontitis. The orange arrows show the area of the periapical lesion in the sagittal (**3A** and **3C**) and coronal (**3B** and **3D**) planes, visible as a hyperintense area, highlighted when fat-suppression techniques are used (in this case, STIR).

resolution of approximately 0.5 mm, allowing multiplanar reconstruction and integration with digital surgical planning workflows. Fat-suppressed sequences such as STIR allow visualization of inflammatory processes including pulp inflammation and periapical pathology.

In addition to *in vivo* imaging, *ex vivo* experiments using phantoms and human cadaveric specimens have contributed to the understanding of signal characteristics of dental tissues, have helped optimize imaging protocols, and have provided insight into how the sequences could be made better and more efficient [3, 13].

Emerging clinical applications

Several possible clinical applications for ddMRI are currently being explored. Although a well-defined indication for the cases in which ddMRI would be a “must have” is currently lacking, research is progressing quickly on defining the cases to which ddMRI would add considerable diagnostic value.

In endodontics, ddMRI may allow assessment of pulp vitality and visualization of periapical inflammation [9]. In orthodontics, ddMRI has been investigated as a potential radiation-free alternative for cephalometric analysis. A pilot study demonstrated high reliability for identifying craniofacial landmarks on ddMRI datasets [14]. In oral surgery, ddMRI may improve visualization of anatomical structures such as the mandibular canal and surrounding soft tissues, potentially improving surgical planning for procedures such as third molar removal and implant placement [10]. MRI is already the reference standard for temporomandibular joint soft-tissue imaging, and ddMRI may allow simultaneous evaluation of TMJ structures and dental anatomy within a single examination [10].

Ongoing studies are addressing other possible uses of this novel imaging modality, specifically in the areas of dental implantology, oral rehabilitation, detection of osteonecrosis of the jaw, and trauma (especially in pediatric patients). In the coming years, results from these studies will help clarify the added value that ddMRI can bring to these areas of dentistry.

Future challenges and research priorities

Despite promising early results, a number of challenges must be addressed to increase the added diagnostic value of this novel modality, and to enable ddMRI to become widely integrated into routine dental practice.

First, clinicians require training to interpret ddMRI images of dentomaxillofacial structures, as these differ substantially from conventional radiographic images (e.g., panoramic images and CBCT volumes).

Second, even though image acquisition times are short compared to “regular” MRI, they are still longer than those of conventional, ionizing-radiation-based imaging modalities (e.g., minutes instead of seconds). Emerging

technologies, including compressed sensing and AI-based reconstruction, may reduce scan times in future system upgrades, making them an even better fit for dental settings, where several patients must be imaged per hour.

Third, susceptibility artifacts from metallic restorations are a technical challenge that needs to be addressed from multiple perspectives [6]. While the lower field strength causes artifacts to be milder than those usually present in “regular” MRI, some dental materials (e.g., stainless-steel wires from orthodontic appliances) still compromise diagnostic image quality in their vicinity. Within this context, the causal relationship between diverse dental materials, ddMRI artifacts, and deterioration in diagnostic image quality must be investigated, and the technical possibilities for mitigating these artifacts must be further explored.

Finally, multicenter clinical validation studies will be required to establish diagnostic accuracy, clinical indications, and the cost-effectiveness of ddMRI compared with existing imaging modalities. We understand that this step will probably become routine in the coming years, when other schools of dentistry around the globe join the ddMRI initiative to reduce the use of ionizing radiation while adding diagnostic value to the images that dentists regularly acquire.

Conclusion

The early experience with ddMRI in the Department of Dentistry and Oral Health at Aarhus University has shown that this technology has the potential to expand diagnostic imaging capabilities in dentistry. By enabling radiation-free imaging with simultaneous visualization of hard and soft tissues, ddMRI is a promising complement to conventional radiographic techniques. Continued collaboration between dentistry, radiology, and MRI physics, as well as close partnerships between the life-science industry and universities will be essential to further develop this modality and determine its role in future dentomaxillofacial imaging.

References

- 1 Tyndall DA, Price JB, Gaalaas L, Spin-Neto R. Surveying the landscape of diagnostic imaging in dentistry's future: Four emerging technologies with promise. *J Am Dent Assoc.* 2024;155(5):364–378. doi: 10.1016/j.adaj.2024.01.005.
- 2 Fuglsig JMCES, Reis INRD, Yeung AWK, Bornstein MM, Spin-Neto R. The current role and future potential of digital diagnostic imaging in implant dentistry: A scoping review. *Clin Oral Implants Res.* 2024;35(8):793–809. doi: 10.1111/clr.14212.

- 3 Greiser A, Christensen J, Fuglsig JMCES, Johannsen KM, Nixdorf DR, Burzan K, et al. Dental-dedicated MRI, a novel approach for dentomaxillofacial diagnostic imaging: technical specifications and feasibility. *Dentomaxillofac Radiol.* 2024;53(1):74–85. doi: 10.1093/dmfr/twad004.
- 4 Hovener JB, Zwick S, Leupold J, Eisenbeibeta AK, Scheifele C, Schellenberger F, et al. Dental MRI: imaging of soft and solid components without ionizing radiation. *J Magn Reson Imaging.* 2012;36(4):841–6. doi: 10.1002/jmri.23712.
- 5 Demirturk Kocasarac H, Geha H, Gaalaas LR, Nixdorf DR. MRI for Dental Applications. *Dent Clin North Am.* 2018;62(3):467–480. doi: 10.1016/j.cden.2018.03.006.
- 6 Johannsen KM, Christensen J, Matzen LH, Hansen B, Spin-Neto R. Interference of titanium and zirconia implants on dental-dedicated MR image quality: *ex vivo* and *in vivo* assessment. *Dentomaxillofac Radiol.* 2025;54(2):132–139. doi: 10.1093/dmfr/twae071.
- 7 Flugge T, Gross C, Ludwig U, Schmitz J, Nahles S, Heiland M, et al. Dental MRI-only a future vision or standard of care? A literature review on current indications and applications of MRI in dentistry. *Dentomaxillofac Radiol.* 2023;52(4):20220333. doi: 10.1259/dmfr.20220333.
- 8 Liedke GS, Christensen J, Spin-Neto R. Evaluation of Temporomandibular Joint Structures and Conditions Using Dental-Dedicated MRI: A Feasibility Study. *J Oral Rehabil.* 2026;53(1):149–159. doi: 10.1111/joor.70071.
- 9 Johannsen K, Christensen J, Kirkevang LL, Matzen LH, Spin-Neto R. Dental-Dedicated MRI in the Assessment of Pulp Vitality and Apical Periodontitis: A Feasibility Study. *Int Endod J.* 2026. doi: 10.1111/iej.70101. Epub 2026 Jan 13.
- 10 Christensen J, Spin-Neto R, Matzen LH. Comparison of dental-dedicated MRI to 2D radiographic images and cone beam CT in the assessment of lower third molars: a prospective study. *Dentomaxillofac Radiol.* 2026. doi: 10.1093/dmfr/twag011. Epub ahead of print.
- 11 Nixdorf DR, Greiser A, Hayes C, Gaalaas L, Groenke BR, Fuglsig JMCES, et al. Comparison of a 0.55 T dental-dedicated magnetic resonance imaging system with a 1.5 T system in evaluation of the temporomandibular joint regarding subjective image quality assessment and rater agreement. *Oral Surg Oral Med Oral Pathol Oral Radiol.* 2025;140(1):113–124. doi: 10.1016/j.oooo.2025.02.011. Epub 2025 Feb 19.
- 12 Fuglsig JMCES, Sampaio-Oliveira M, Spin-Neto R. Dental-dedicated magnetic resonance imaging in the follow-up of lower third molar removal. *Oral Radiol.* 2025;41(2):296–301. doi: 10.1007/s11282-024-00787-x.
- 13 Sampaio-Oliveira M, Oliveira ML, Spin-Neto R. Development of a carrier device for dental-dedicated magnetic resonance imaging. *Imaging Sci Dent.* 2025;55(2):207–13. doi: 10.5624/isd.20240254.
- 14 Srivastav S, Spin-Neto R, Aiyar A, Stoustrup PB. Accuracy and reliability of magnetic resonance imaging in orthodontic diagnosis and treatment planning-a systematic review and meta-analysis. *Eur J Orthod.* 2024;46(3):cjae019. doi: 10.1093/ejo/cjae019.

Contact

Rubens Spin-Neto, Dr. odont., Ph.D. in Oral Radiology, Ph.D. in Periodontology, M.S. Professor in Dentomaxillofacial Radiology and Imaging
Department of Dentistry and Oral Health
Aarhus University
Vennelyst Boulevard 9
8000 Aarhus
Denmark
Tel.: +45 42 43 48 36
rsn@dent.au.dk



MRI for the Dental Practice. A Technology Finds its Way

Miriam Keil

Siemens Healthineers, Erlangen, Germany

Magnetic resonance imaging (MRI) is an indispensable tool in the field of radiology. For decades, our MRI systems have been supporting diagnosis and treatment in clinical fields such as orthopedics, neurology, and oncology, to name just a few.

To date, MRI has played only a marginal role in dentistry. Yet as a non-ionizing imaging modality that provides superior soft tissue contrast, MRI has much to offer this discipline. Onsite access to an MRI system would help dentists fill the diagnostic gap left by their traditional dental imaging systems.

When we started exploring how we could bring MRI into dental practices, it quickly became clear that a dental-dedicated MRI (ddMRI) system would have to fulfil some very specific requirements: From our ethnographic research, we gained a clear picture of the unique operational and infrastructural realities in dental practices. These insights led us to fundamentally rethink MRI and create a system that is tailored precisely to the workflows, personnel, and clinical questions in dentistry. The aim was to enable dentists, their teams, and their patients to benefit from imaging that shows, for instance, whether a nerve is too close to a wisdom tooth, whether pulp is vital, or whether the temporomandibular joint (TMJ) is structurally damaged.

The diagnostic gap


To understand why MRI matters for dentistry, it helps to look at what radiologists have been using routinely for decades. In orthopedics, no cruciate ligament tear is diagnosed without MRI; no herniated disc, no meniscal lesion, no rotator cuff rupture. In neurology, MRI is the gold standard for virtually every question. In oncology, cardiology, gastroenterology — wherever soft tissue contrast is decisive — MRI has become indispensable.

Dentistry, by contrast, still works primarily with ionizing radiation. The orthopantomogram shows bone and teeth in an overview image. Cone beam computed tomography (CBCT) delivers three-dimensional representations with excellent resolution for hard tissues. Both of these modalities are invaluable and firmly established in dental practice. But they share a fundamental limitation: Soft tissue appears largely uniform, with little structural differentiation.


In practical terms, this means that the TMJ articular disc (a cartilaginous structure only a few millimeters thick and a frequent cause of TMJ disorders) is invisible on CBCT. The dental pulp, whose vitality may determine whether a tooth can be preserved, can only be assessed indirectly. The inferior alveolar nerve — which must be protected during implant placement or wisdom tooth extraction — can only be inferred from the surrounding bony canal, not visualized directly.

This diagnostic gap is no secret. Dentists know there are clinical questions for which their available imaging modalities provide no satisfactory answers. Until recently, the solution has been to refer patients to radiology practices or hospitals for MRI examinations in special cases. This is cumbersome, time-consuming, and simply not practical for many clinical questions.


Tooth extraction workflow

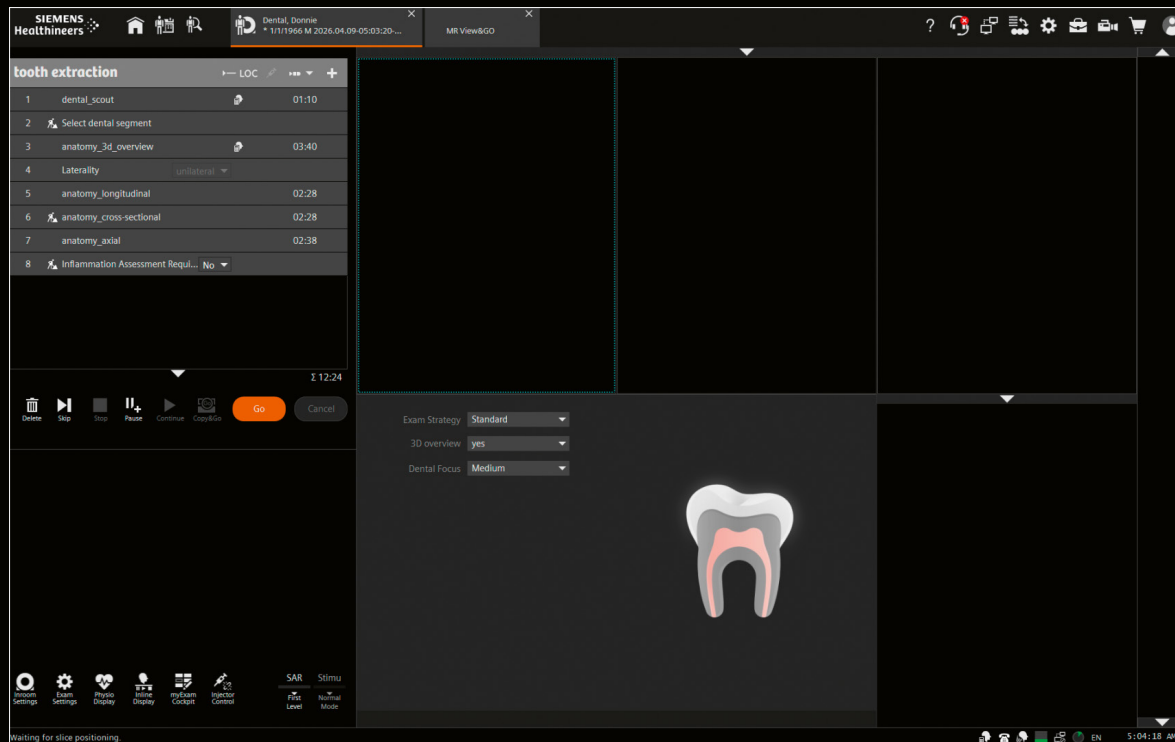

- 1 Tooth extraction is the most complex workflow. It requires three primary user interventions (indicated by the  operator icon): dental segment selection, fine adjustment of cross-sectional slices, and optional inclusion of additional inflammation sequences.

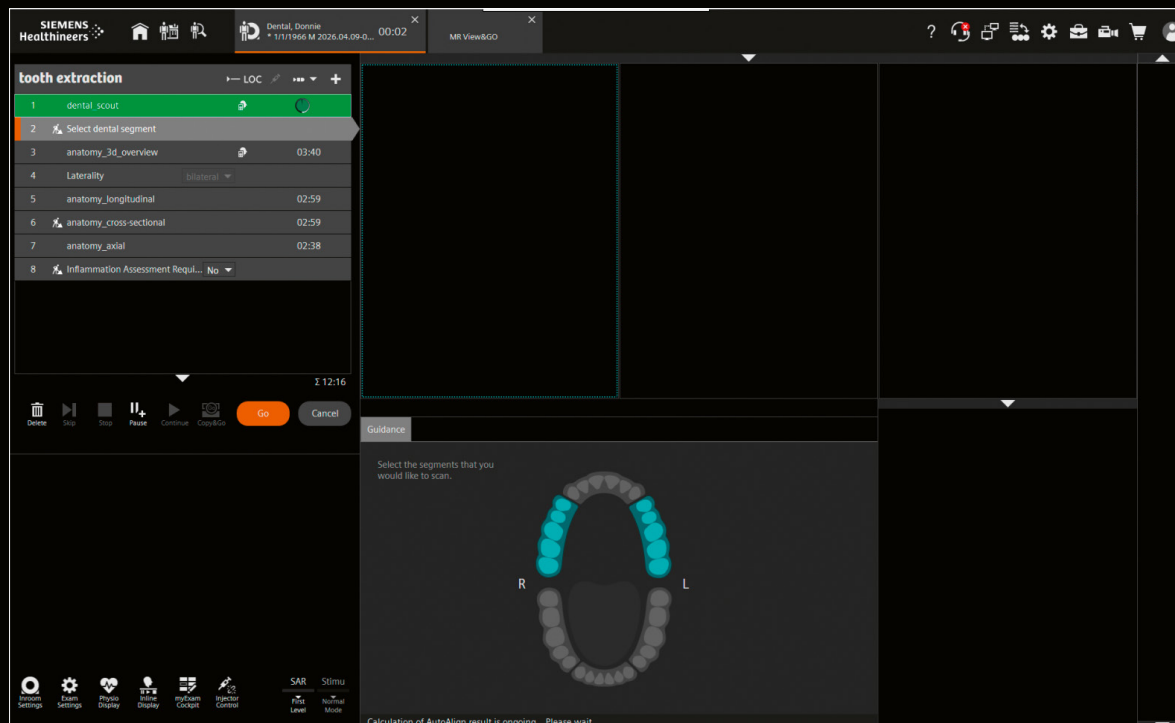

- 2 The selection of dental segments relevant for the examination follows established cone beam CT workflows. This supports intuitive use and easy recognition for users familiar with CBCT.

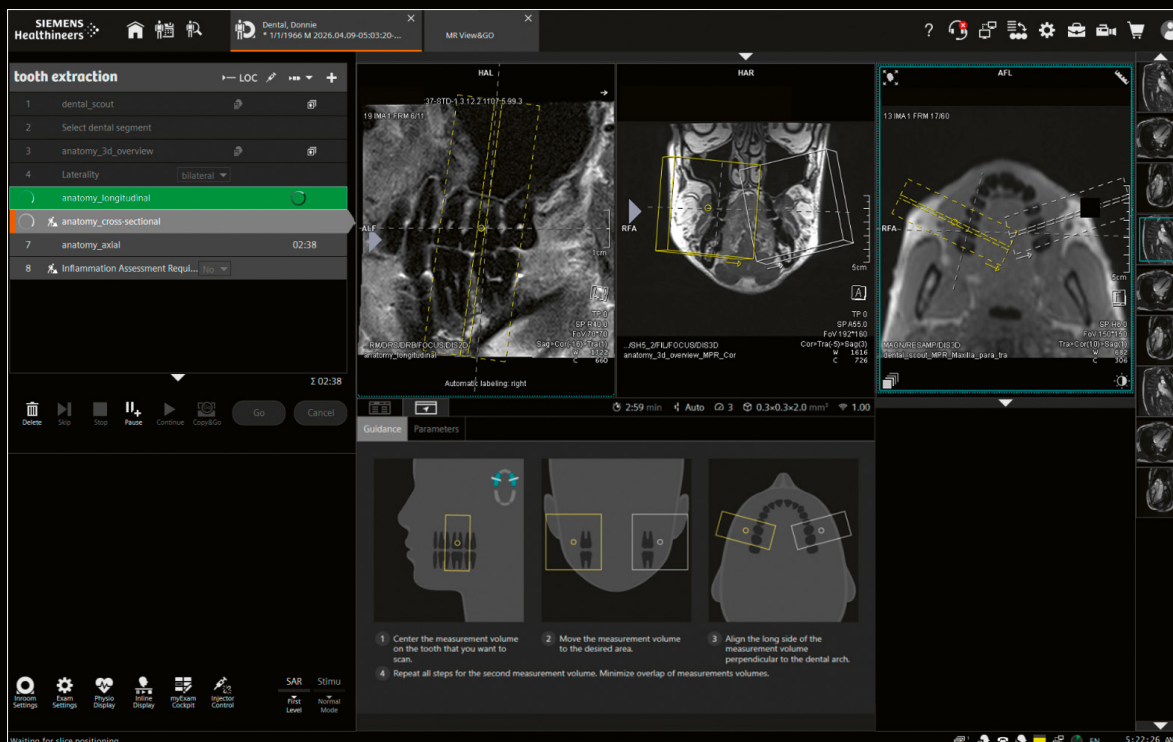


- 3 The only interactive planning step occurs on the cross-sectional images. This is to achieve correct tooth selection and adequate coverage of the root or area of interest while minimizing acquisition time. The user is guided through this process.



- 4 If inflammation is suspected, the user can extend the workflow with additional steps, which are automatically planned in advance.

TMJ workflow

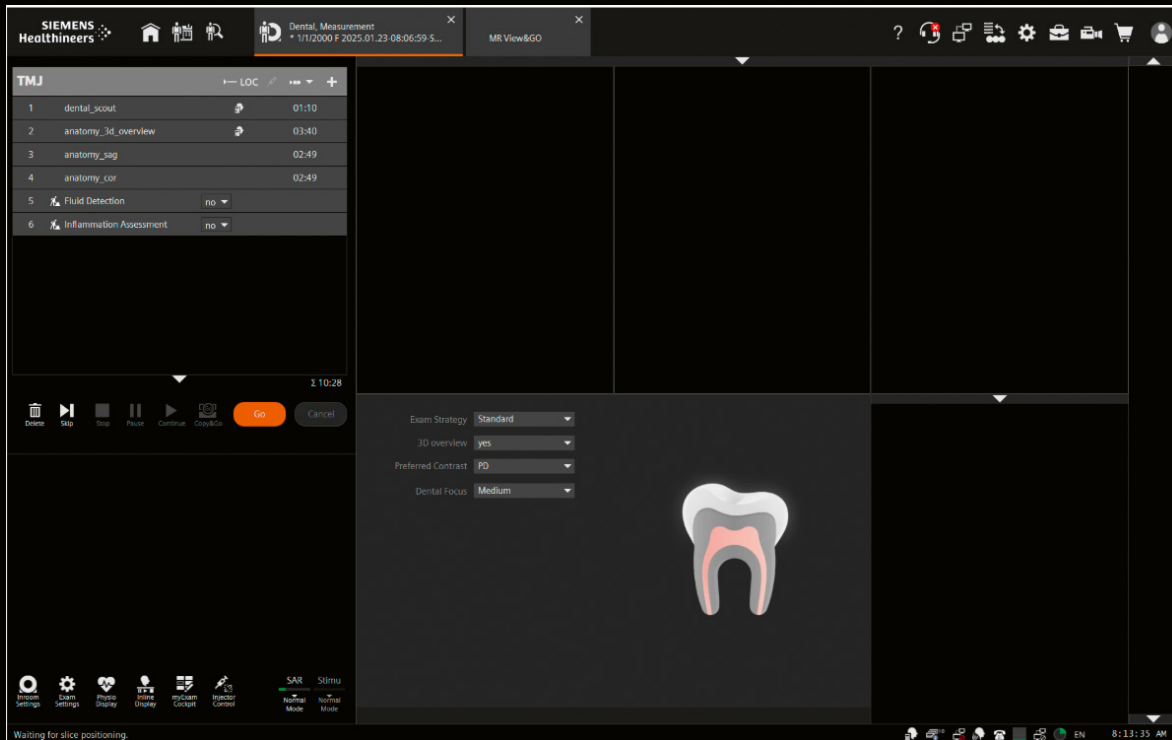

- 5 The TMJ workflow is fully automated. User interaction is only required if additional contrast acquisitions for fluid detection or inflammation assessment are required.

Why a conventional MRI scanner is not the answer

At first glance, the solution might seem obvious: Install a conventional MRI system in the dental practice. In reality, this won't work due to fundamental workflow differences. Conventional MRI systems are designed for operation by specialist personnel, who make decisions about anatomical planes, sequence types, slice thicknesses, repetition and echo times, flip angles, and fat suppression. Each parameter influences image quality and diagnostic value. Radiology technologists spend years becoming experts in these relationships.

Dental practices do not have this workforce — and maintaining one would make no economic sense. A dental assistant who provides chairside assistance, schedules appointments, and takes X-rays cannot be expected to also become an MRI specialist on the side. Dental practices therefore need an MRI system that non-specialist staff can operate without compromising image quality or outcomes.

Spatial and infrastructural constraints are another issue. Conventional MRI systems require large rooms, elaborate shielding, a helium supply with quench pipes, and reinforced floors to take the magnet's weight. Most dental practices lack both the space and the financial means to make such modifications.

From observation to solution

The development of our ddMRI system — MAGNETOM Free.Max Dental Edition — began not with technical specifications, but with ethnographic field research. We shadowed dentists, analyzed workflows, and spoke extensively with practitioners. We learned that a 45-minute MRI exam is incompatible with 15-minute appointment slots; that complex planning steps will not survive daily practice; that dentists think in tooth regions and jaw segments, not coronal and sagittal planes; that they care exclusively about the dentomaxillofacial region; and that technical terminology creates unnecessary barriers.

From these insights, clear design principles emerged: Examinations should be indication-driven, not parameter-driven. Instead of selecting sequences, users choose clinical questions or anatomical targets. For each indication — whether TMJ diagnostics, nerve visualization, pulp assessment, or periodontal evaluation — the system automatically configures all technical settings. The myExam Dental Assist software guides the user through every step, from patient positioning to the finished scan. In routine cases, the full examination can be completed in less than 15 minutes.

Patient positioning is handled via Select&GO with automatic isocentering. This is performed at the touch of a button instead of via manual adjustments. The Dental Segment Planning Add-In extends automation to tooth-specific imaging: Clinicians simply select the relevant dental segment, then the system identifies the anatomy, plans slice groups, and sets imaging geometry automatically using AutoAlign. Cross-sectional planes can be fine-tuned for precise local imaging, with support from dedicated guidance tools. The system enables TMJ examinations to be planned fully automatically.

Imaging quality is supported by Deep Resolve, an AI-based reconstruction technology that achieves 0.3-millimeter in-plane resolution with substantially reduced scan times. Examination programs carry descriptive dental names rather than technical sequence labels, bridging the language gap between MRI technology and dental practice.

We also had to rethink the hardware. The standard head coils used in radiology are optimized for broad anatomical coverage — from the vertex to the base of the neck. For dental imaging, they are too large, difficult to position, and unsuitable for dedicated dentomaxillofacial imaging. The newly developed 7-channel dental coil features an open architecture for patient comfort and access. It can be secured with one hand and is adjustable in the tilt angle and longitudinal direction. Reproducible alignment is facilitated by integrated positioning markers at the glabella and orbits. The dental coil also enables the use of established dental positioning aids, such as bite wedges. The scanner's 80-centimetre bore diameter, open coil design, and dedicated cushioning system substantially improve patient comfort, which is particularly relevant for claustrophobic patients.

Seeing more

Achieving clinically useful image quality within dental time constraints relies on advanced reconstruction technology. Deep Resolve is a deep-learning-based algorithm trained on high-quality MRI datasets to recover spatial resolution and reduce noise from undersampled acquisitions.

Traditional MRI balances resolution, signal-to-noise ratio, and acquisition time. High-resolution imaging usually means long scans or noisy images. Deep Resolve relaxes this trade-off by reconstructing high-quality images from less data. The result is up to 60% shorter scan times, while maintaining 0.3-millimeter in-plane resolution.

This resolution is enough for visualizing structures that are critical to dental diagnosis: the pulp chambers, root canals, alveolar cortex, and the TMJ disc.

While CBCT still surpasses MRI for pure bone morphology, ddMRI offers a different diagnostic value: Its strength lies in its ability to visualize soft tissue contrast, inflammation, vascularization, and pulp vitality. This is often more decisive than incremental gains in geometric resolution.

Optimized soft tissue contrast enables clear differentiation of the pulp-dentin interface, inferior alveolar and lingual nerves, and soft tissue components of the TMJ. These capabilities fundamentally exceed what radiographic techniques can provide, and they complement CBCT rather than replace it.

Focusing on what matters

A unique challenge in dental MRI is clinical focus. Dentists are experts in dentomaxillofacial anatomy. Intracranial structures are usually irrelevant and introduce distraction rather than clarity.

The Dental Focus addresses this through acquisition-level field-of-view restriction. Three calibrated modes (narrow, medium, wide) ensure that only clinically relevant anatomy is displayed at diagnostic quality. The result is focused, unambiguous imaging tailored to dental decision-making.

We also adapted the orientations and the nomenclature of MRI. Instead of coronal and sagittal views, images are presented using dental orientations: longitudinal, cross sectional, and axial views aligned to the teeth and jaws. AutoAlign algorithms based on cephalometric landmarks provide consistent orientation. The system is capable of producing three-dimensional panoramic reconstructions, longitudinal sections aligned with the dental arch, and cross-sectional views perpendicular to specific teeth.

An MRI system that fits into dental practices

The MAGNETOM Free.Max platform removes infrastructural barriers. It operates with only 0.7 liters of helium, which is sealed in a maintenance-free system that requires no refilling and no quench pipe. Space requirements are under 25 square meters. At less than 3.2 tons and under two meters in height, the system also fits through standard doors. Siting complexity is therefore dramatically reduced compared to high-field MRI systems.

The diagnostic scope of dental-dedicated MRI

Clinical applications span multiple domains of oral and maxillofacial medicine. In TMJ disorders, MRI is the only non-invasive modality that directly visualizes the articular disc, showing its position, shape, and integrity. Differentiating between muscular dysfunction, disc displacement with or without reduction, and degenerative joint disease is clinically decisive and fundamentally influences the therapeutic approach.

In endodontics, MRI enables direct visualization of the pulp. Current vitality tests — thermal or electrical — provide indirect, sometimes ambiguous information. MRI can display inflammatory changes around the apex, enable assessments of vascularization, and support the detection of periapical pathology at early stages before substantial bone loss has occurred.

Periodontics benefits from visualization of soft tissue inflammation and ligament integrity. Active periodontitis shows characteristic edema and hyperemia on MRI before radiographically detectable bone loss occurs. The extent of soft tissue attachment loss can be assessed non-invasively, complementing clinical probing.

MRI makes surgical planning more precise, thanks to the detailed nerve visualization. The inferior alveolar nerve, mental nerve, and lingual nerve can be traced along their courses. This allows more precise planning of implant placement, wisdom tooth extractions, and orthognathic procedures, particularly in cases where the nerve course is unusually close to the planned surgical site or where previous interventions have altered normal anatomy.

In juvenile dentistry, the absence of ionizing radiation is particularly relevant. Repeated MRI in growing children — to monitor dental development, assess impacted teeth, or follow trauma over time — carries no cumulative radiation burden.

Beyond these applications, the technology opens up possibilities for clinical questions that have traditionally been difficult to address. Osseous lesions with predominantly medullary involvement can be detected earlier by MRI than by CT. Soft tissue tumors and cystic lesions can be differentiated more confidently. Osteonecrosis of the jaw — whether medication-related or from other causes — shows characteristic marrow signal changes on MRI before radiographically visible bone destruction occurs.

A system tailored to the practical realities of dentistry

With ddMRI, magnetic resonance imaging enters dentistry not as a tool simply transplanted from radiology, but as a system designed for the realities of dentistry.

What makes this transition viable is not mere miniaturization or cost reduction, but rather comprehensive adaptation. The ddMRI system works because it respects the constraints of dental practice: limited time per patient, operation by personnel without specialist training, integration with existing imaging modalities, and strict adherence to scope-of-practice boundaries. It delivers genuine diagnostic value — soft tissue visualization that cannot be obtained by other means — while fitting into established workflows.

Dentists gain access to a diagnostic modality that was previously out of reach. Patients benefit from improved diagnostics without referrals or radiation exposure. The technology makes visible what has long been invisible.



Contact

Dr. rer. nat. Miriam Keil
Siemens Healthineers
DI MR NEW
Doris-Ruppenstein-Str. 4
91052 Erlangen
Germany
miriam.keil@siemens-healthineers.com

Not for distribution in the US

On account of certain regional limitations of sales rights and service availability, we cannot guarantee that all products included in this brochure are available through the Siemens Healthineers sales organization worldwide. Availability and packaging may vary by country and is subject to change without prior notice. Some/All of the features and products described herein may not be available in the USA.

The information in this document contains general technical descriptions of specifications and options as well as standard and optional features which do not always have to be present in individual cases, and which

may not be commercially available in all countries. Due to regulatory reasons their future availability cannot be guaranteed. Please contact your local Siemens Healthineers organization for further details.

Siemens Healthineers reserves the right to modify the design, packaging, specifications, and options described herein without prior notice. Please contact your local Siemens Healthineers sales representative for the most current information.

Note: Any technical data contained in this document may vary within defined tolerances. Original images always lose a certain amount of detail when reproduced.

Siemens Healthineers Headquarters

Siemens Healthineers AG
Siemensstr. 3
91301 Forchheim, Germany
Phone: +49 9191 180
siemens-healthineers.com